

LOAN DOCUMENT

PHOTOGRAPH THIS SHEET

DTIC ACCESSION NUMBER

LEVEL

INVENTORY

197 - 1995

DOCUMENT IDENTIFICATION

MAR 1995

DISTRIBUTION STATEMENT A

Approved for public release;
Distribution Unlimited

DISTRIBUTION STATEMENT

ACCESSION FOR

NTIS GRA&I
DTIC TRAC
UNANNOUNCED
JUSTIFICATION

DDI

BY

DISTRIBUTION/

AVAILABILITY CODES

DISTRIBUTION

AVAILABILITY AND/OR SPECIAL

A-1

DISTRIBUTION STAMP

19980731 120

DATE RECEIVED IN DTIC

REGISTERED OR CERTIFIED NUMBER


PHOTOGRAPH THIS SHEET AND RETURN TO DTIC-FDAC

H
A
N
D
L
E

W
I
T
H

C
A
R
E

THE QUANTITATIVE ION EXCHANGE SEPARATION OF URANIUM FROM IMPURITIES

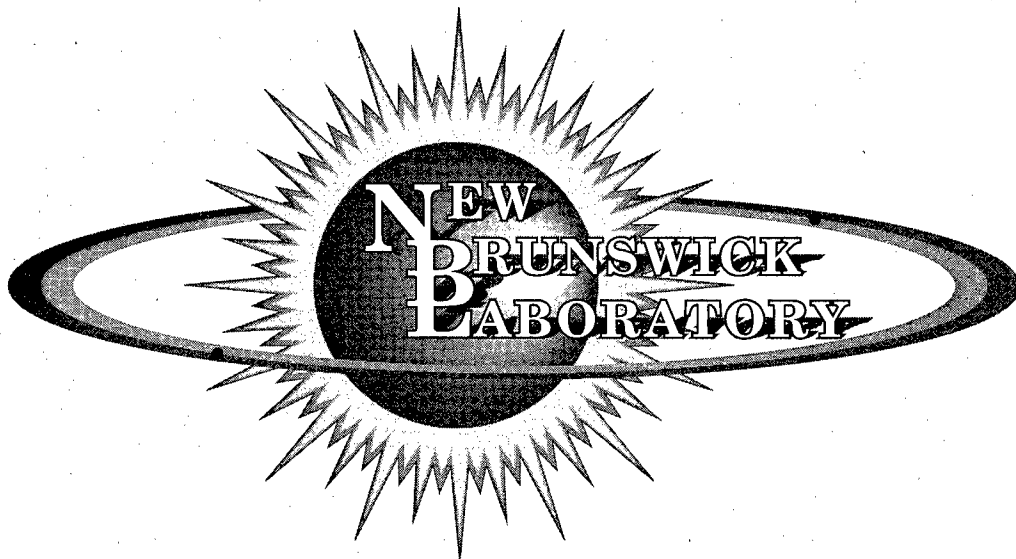

**Usha I. Narayanan, Peter B. Mason, Jeffrey P. Zebrowski,
Matthew Rocca, Iris W. Frank, Marianne M. Smith,
Kimberly D. Johnson, Glennda J. Orlowicz
and Eric Dallmann**



REPORT DOCUMENTATION PAGEForm Approved
OMB No. 074-0188

Public reporting burden for this collection of information is estimated to average 1 hour per response, including the time for reviewing instructions, searching existing data sources, gathering and maintaining the data needed, and completing and reviewing this collection of information. Send comments regarding this burden estimate or any other aspect of this collection of information, including suggestions for reducing this burden to Washington Headquarters Services, Directorate for Information Operations and Reports, 1215 Jefferson Davis Highway, Suite 1204, Arlington, VA 22202-4302, and to the Office of Management and Budget, Paperwork Reduction Project (0704-0188), Washington, DC 20503

1. AGENCY USE ONLY (Leave blank)		2. REPORT DATE March 1995	3. REPORT TYPE AND DATES COVERED Research and Development Report, March 1995	
4. TITLE AND SUBTITLE The Quantitative Ion Exchange Separation of Uranium From Impurities			5. FUNDING NUMBERS N/A	
6. AUTHOR(S) Usha I. Narayanan, Peter B. Mason, Jeffrey P. Zebrowski, Matthew Rocca, Iris W. Frank, Marianne M. Smith, Kimberly D. Johnson, Glenna J. Orlowcz and Eric Dallmann				
7. PERFORMING ORGANIZATION NAME(S) AND ADDRESS(ES) New Brunswick Laboratory Argonne, Illinois			8. PERFORMING ORGANIZATION REPORT NUMBER NBL-331	
9. SPONSORING / MONITORING AGENCY NAME(S) AND ADDRESS(ES) SERDP 901 North Stuart St. Suite 303 Arlington, VA 22203			10. SPONSORING / MONITORING AGENCY REPORT NUMBER N/A	
11. SUPPLEMENTARY NOTES Prepared for the U.S. Department of Energy. This work was supported in part by SERDP. The United States Government has a royalty-free license throughout the world in all copyrightable material contained herein. All other rights are reserved by the copyright owner.				
12a. DISTRIBUTION / AVAILABILITY STATEMENT Approved for public release: distribution is unlimited			12b. DISTRIBUTION CODE A	
13. ABSTRACT (Maximum 200 Words) Two methods were tested for the quantitative separation of uranium from elemental impurities using commercially available resins. The sorption and elution behavior of uranium and the separation of it from a variety of other elements was studied. The first method utilized an anion exchange resin while the second method employed an extraction resin. The first method, the anion exchange of uranium (VI) in an acid chloride medium, was optimized and statistically tested for quantitative recovery of uranium. This procedure involve adsorption of uranium onto Bio-Rad AG 1X8 or MP-1 ion exchange resins in 8 M HCl, separation of uncomplexed or weakly complexed matrix ions with an 8 M HCl wash, and subsequent elution of uranium with 1 M HCl. Matrix ions more strongly adsorbed than uranium were left on the resin. Uranium recoveries with this procedure averaged greater than 99.9% with a standard deviation of 0.1%. In the second method, recovery of uranium on the extraction resin did not meet the criteria of this study and further examination was terminated.				
14. SUBJECT TERMS uranium, extraction resin, Bio-Rad AG 1-X8, MP-1, SERDP			15. NUMBER OF PAGES 66	
			16. PRICE CODE N/A	
17. SECURITY CLASSIFICATION OF REPORT unclass	18. SECURITY CLASSIFICATION OF THIS PAGE unclass	19. SECURITY CLASSIFICATION OF ABSTRACT unclass	20. LIMITATION OF ABSTRACT UL	

NSN 7540-01-280-5500

Standard Form 298 (Rev. 2-89)
Prescribed by ANSI Std. Z39-18
298-102

DTIC QUALITY INSPECTED 1

U.S. DEPARTMENT OF ENERGY

THE QUANTITATIVE ION EXCHANGE SEPARATION OF URANIUM FROM IMPURITIES

Usha I. Narayanan, Peter B. Mason, Jeffrey P. Zebrowski,
Matthew Rocca, Iris W. Frank, Marianne M. Smith,
Kimberly D. Johnson, Glenda J. Orlowicz and Eric Dallmann

Research and Development Report

NEW BRUNSWICK LABORATORY

Argonne, Illinois

March, 1995

H. Rodney Martin, Acting Director

NOTICE

REFERENCE TO A PARTICULAR BRAND OR MANUFACTURER OF EQUIPMENT DOES NOT CONSTITUTE AN ENDORSEMENT BY THE NEW BRUNSWICK LABORATORY OR THE DEPARTMENT OF ENERGY.

TABLE OF CONTENTS

ACKNOWLEDGMENTS	v
ABSTRACT	vii
INTRODUCTION	1
EXPERIMENTAL	2
Apparatus, Materials and Reagents	2
PROCEDURE	2
Test Samples	2
Column Capacities	3
Column Preparation	3
Sample Preparation	4
Sample Loading	4
Column Wash	4
Uranium Elution	5
Uranium and Impurity Analysis	5
RESULTS AND DISCUSSION	6
Sorption and Elution Profiles of Pure Uranium Solutions with U/TEVA	7
Table I - Elution Pattern of Uranium on U/TEVA Column	8
Table II - Sorption and Elution Profile of U/TEVA resin with 0.01 M Ammonium Bioxalate Eluant	9
Sorption and Elution Profiles of Pure Uranium Solutions with Bio-Rad Resins	10
Recovery of Pure Uranium with Bio-Rad Resins	10
Table III - Adsorption and Elution Behavior of Uranium on Bio-Rad Resins	11
Table IV - Uranium Recovery with Bio-Rad Resins	12
Table V - Optimized Uranium Recovery with Bio-Rad 1-X8 and MP-1 Resins	13
Ion Exchange of the Matrix Blank	13
Ion-Exchange of Spiked Uranium Samples	14
Table VI - Sorption and Elution Pattern of 22 Impurities on Bio-Rad MP-1 Resin	15
Table VII - Uranium Recovery of Spiked Samples	17
Uranium Breakthrough and Recovery	18
Statistical Testing	18

Table VIII - Uranium Breakthrough Study on Bio-Rad AG 1-X8 and MP-1 Resins	19
Table IX - Summary Table of Resin Statistics for Uranium Recovery	20
CONCLUSIONS	21
REFERENCES	21
APPENDIX A: Procedure for the Purification of Uranium by Anion Exchange	23
APPENDIX B: Statistical Plan of Analysis	37
APPENDIX C: Statistical Analysis of Uranium Purification Project	43
DISTRIBUTION LIST	55

ACKNOWLEDGMENTS

This work was performed under the auspices of the Office of Safeguards and Security in the Office of Non-Proliferation and National Security. Funding for this work was provided by the Strategic Environmental Research and Development Program (SERDP) under the Pollution Prevention thrust area.

The authors express their appreciation to R. D. Oldham for providing technical guidance, to M. I. Spaletto for technical discussions, to P. V. Croatto for providing assistance during KPA analyses, to D. T. Baran, G. A. Sowell, L. D. Preston and A. V. Stiffin for X-ray analyses, to W. G. Mitchell and M. I. Spaletto for reviewing this report, and to M. Ficner for secretarial support.

ABSTRACT

Two methods were tested for the quantitative separation of uranium from elemental impurities using commercially available resins. The sorption and elution behavior of uranium and the separation of it from a variety of other elements was studied. The first method utilized an anion exchange resin while the second method employed an extraction resin. The first method, the anion exchange of uranium (VI) in an acid chloride medium, was optimized and statistically tested for quantitative recovery of uranium. This procedure involved adsorption of uranium onto Bio-Rad AG 1-X8 or MP-1 ion exchange resins in 8 *M* HCl, separation of uncomplexed or weakly complexed matrix ions with an 8 *M* HCl wash, and subsequent elution of uranium with 1 *M* HCl. Matrix ions more strongly adsorbed than uranium were left on the resin. Uranium recoveries with this procedure averaged greater than 99.9% with a standard deviation of 0.1%. In the second method, recovery of uranium on the extraction resin did not meet the criteria of this study and further examination was terminated.

INTRODUCTION

The investigation into a separation procedure with quantitative recovery of uranium was prompted by a need to meet waste reduction requirements and reference material characterization demands. At New Brunswick Laboratory (NBL), a major emphasis has been placed on developing uranium assay measurement methods which will avoid the generation of mixed wastes. The measurement methods under development at NBL, however, are not expected to be as free of matrix interferences as existing procedures and therefore will require prior removal of impurities. In addition, the necessity for a high degree of accuracy in the uranium assays performed at NBL demands a purification procedure providing essentially quantitative (>99.9%) recovery of uranium. Therefore, uranium separation procedures that provide large separation factors for a variety of impurities were selected and investigated for quantitative recovery of uranium.

A literature survey of a variety of uranium separation methods was performed.^{1,2} Based on past experience and success with developing a plutonium ion exchange method at NBL,^{3,4} procedures based on the adsorption of anionic⁵ and cationic⁶ complexes of uranium and sorption by extraction resins⁷ were reviewed. The separation methods selected for evaluation had to have the potential to separate 27 specific impurities expected to be present in one or more of the variety of sample types received at NBL.⁸ Another consideration for selection was the type and quantity of waste generated in the separation procedure. From this survey, two methods were selected for laboratory testing.

The first method involved the anion exchange of U(VI) species in chloride media. In the chosen procedure, uranium was sorbed as the chloride complex in 8 M HCl onto strongly basic anion exchange resin, separated from unretained elements by washing the resin with 8 M HCl, then selectively eluted in 1 M HCl. The second method employed U/TEVA resin, a resin coated with a complexing agent that has a very high specificity for uranium. In this method, the sample was loaded in 2.5 M HNO₃, impurities were separated with a 2.5 M HNO₃ wash, and uranium was eluted in 0.01 M HNO₃. The U/TEVA method had been successfully tested at NBL and used for the purification of uranium in preparation for mass spectrometric analysis.⁹

The variables affecting the sorption and elution of the various sorbed complexes and the recovery of uranium, including acid concentrations of the load and wash solutions and the degree of cross-linking of the ion exchange resin, were investigated and are described in this report. Statistical testing of an optimized procedure using Bio-Rad AG 1-X8 and AG MP-1 resins is also discussed. The sample size tested for quantitative recovery of uranium with the anion exchange method was 12 mg. The effect of impurities on the quantitative recovery of 12 mg uranium samples was tested over the range of 10% by weight of impurities (mole ratio of U to impurities of 2:1) to 40% by weight of impurities (mole ratio of U to impurities of approximately 1:1.7).

The two methods tested could potentially be applied to a broad spectrum of analytical applications such as the separation of uranium for subsequent determination by constant current coulometry, isotope dilution mass spectrometry, alpha spectroscopy and the removal of possible uranium interference prior to the quantitative determination of impurities by ICP-AES and X-ray fluorescence.

EXPERIMENTAL

Apparatus, Materials and Reagents

Resins: Bio-Rad AG® 1-X2, 1-X4, 1-X8 and MP-1, chloride form, 100-200 mesh (dry) from Bio-Rad Laboratories (Richmond, CA). U/TEVA Spec,® 100-150 μ , Eichrom Industries, Inc. (Darien, IL).

Columns: Bodman Polypropylene Columns, Large Capacity Robotic Tube, Order No. 3300-25; polyethylene frits, 70 μ , Order No. 3201-00, from Bodman Industries (Aston, PA). Total volume is nominally 18 mL; volume of stem that holds the resin is approximately 2.5 mL.

Column Holder: NBL design.

Sample Beaker: Glass and/or PFP Teflon,® 50 mL, 100 mL and 400 mL.

Transfer Pipets: Graduated plastic with an approximate draw of 3.5 mL.

Reagents: Hydrochloric acid, nitric acid and ammonium bioxalate. All chemicals were of ACS grade purity. Deionized distilled water was used for all experiments.

Uranium: NBL Certified Reference Material 112-A, Uranium Metal Assay Standard.

Impurity Spike: Individual standard solutions of 25 elements, all approximately 1000 micrograms of element/mL, atomic spectrometry standard grade, from Inorganic Ventures Inc. (Lakewood, NJ).

PROCEDURE

Test Samples

Uranium samples were prepared by dissolving NBL Certified Reference Material (CRM) 112-A, Uranium Metal Assay Standard, in nitric acid.¹⁰ The subsequent uranyl nitrate solution (2 M HNO₃) was aliquanted into

clean glass and Teflon beakers in appropriate amounts, to deliver the desired 12 mg uranium sample size. Control uranium samples were aliquanted into glass beakers for direct titration by the NBL Titrimetric Method (NBL-Modified Davies and Gray Titration) with no prior ion exchange performed. Elements selected as impurities were combined in definite amounts to create a stock solution of 25 elements. Aliquants of this mixture (consisting of 50 μg , 100 μg or 200 μg each element) were placed into both empty beakers and beakers with uranium to generate matrix blanks and synthetic uranium samples with impurities respectively. All aliquants were dried down on a steam bath and the beakers were covered with parafilm to prevent loss or contamination during storage.

Column Capacities

The column capacities of the U/TEVA and Bio-Rad resins were calculated from experimental values provided by the manufacturers.^{7,11} This information was used to calculate a suitable uranium sample size, small enough to be easily retained on the column, but large enough to be quantitatively analyzed by NBL methods.

The reported capacity of U/TEVA resin is 37 mg U/mL resin, with a recommended working range between 10% and 20% of capacity. Therefore, with a bed volume of 2.5 mL, a sample size of 12 mg uranium represents approximately 15% of resin capacity.

The capacities of the Bio-Rad ion exchange resins were calculated from the exchange capacity data for the resin provided by Bio-Rad in milliequivalent/mL and the charge on the anionic species⁵ in the resin. The capacity of Bio-Rad AG 1-X8 resin is given as 1.4 meq/mL and the anionic species $[\text{H}_3\text{O}(\text{H}_2\text{O})_3(\text{UO}_2\text{Cl}_4)_2]^{3-}$ has a charge of -3. This results in a theoretical capacity of approximately 100 mg U/mL resin. However, for samples that have impurities, the matrix ions also compete for exchange sites. It was decided to work with much lower amounts of uranium and impurities than the calculated column capacities to avoid potential losses. For the desired 2.5 mL resin bed, a 12 mg uranium sample (with no added impurities) represents 5% of the operating capacity of the resin, well below the recommended maximum of 15%.

Column Preparation

The ion exchange column was placed in a column holder above a 100 mL beaker. The resins were used "as is" without any further sieving. A resin slurry was prepared in a clean polyethylene bottle by mixing dry resin with sufficient water to make the slurry. Enough resin was added to the column using a plastic transfer pipet to give a resin height of approximately 4 cm. This resin height provides a wet resin volume of not less than 2.5 mL. To ensure a uniform resin bed, about 5 mL of water was placed in the column, and resin was drawn from the column into the pipet and then released. The resin was allowed to settle in the column.

This helped to remove any air bubbles formed when loading the column, and allowed the resin beads to settle uniformly in the column. The resin bed was conditioned by adding 30 mL of the appropriate acid (8 M HCl reagent for the anion exchange and 2.5 M HNO₃ reagent for the U/TEVA extraction chromatography) to the column with a transfer pipet. For the Bio-Rad AG MP-1 and U/TEVA resins, glass wool was placed at the top of the resin bed to prevent the resin from being disturbed during the addition of reagents. The column beds prepared from other resins remained intact and thus no glass wool was needed for these resins. The flow rate for the U/TEVA column was uniform at 1 mL/min. The average flow rate measured for all of the Bio-Rad resins was 0.6 mL/min. With the AG MP-1 (macroporous) resin, resin fines were observed in the collection beaker during the column preparation step. No visible resin fines were observed in the collection beaker during sample loading or during later steps of the procedure.

Sample Preparation

The Parafilm cover and walls of the sample beaker were rinsed with 12 to 14 mL of the appropriate acid (8 M HCl or 2.5 M HNO₃) and the sample was allowed to dissolve completely.

Sample Loading

The sample solution was slowly and carefully poured, without loss, onto the resin bed. When pouring samples from glass beakers, a small drop of solution tends to adhere to the lip of the beaker. The analyst must be very careful to avoid losing sample at this point. Note that for a 12 mg sample, a loss of just 12 micrograms of uranium represents 0.1% of the sample, thus making recovery of greater than 99.9% impossible. Teflon beakers were also tested to see if their use reduced sample loss during loading since the aqueous solution does not cling to hydrophobic Teflon. A disposable plastic pipet was used to rinse the walls of the sample beaker with 2 mL of acid. The rinse solution was then transferred onto the column. This rinse was repeated four times. The outside of the beaker spout was then rinsed into the column with a few drops of acid. The column effluent from the sample loading was discarded or analyzed for trace uranium or impurities, as appropriate.

Column Wash

Aliquots of 5 mL of acid wash were added to the column and allowed to drain to the top of the resin bed. A total volume of 20 mL wash was utilized in the final anion exchange procedure. The effluent (washings) was discarded or set aside for appropriate analysis (uranium and/or impurities).

Uranium Elution

Uranium was eluted into a clean labeled beaker by adding the first 5 mL of the appropriate eluting solution (0.01 M HNO₃ or 1 M HCl) to the column dropwise. The dropwise additions were used to wash the sides of the reservoir. Eluting solution was then added in aliquots of 5 mL, allowing the column to drain between increments, until a total of 40 mL was added. Small fractions of the eluate were collected and analyzed to determine elution curves. The final volume of eluant for the ion exchange procedure was determined to be 40 mL based upon analysis of the elution curves.

Uranium and Impurity Analysis

The eluate was evaporated on a hot plate and fumed twice with 2 mL of concentrated HNO₃. The collections from sample loadings and washes were also evaporated and fumed as described above. The sample loadings, washes, and eluate fractions were analyzed for uranium and impurities, as appropriate. For these experiments, three analytical methods were used for the assay of uranium and impurities.

A laser kinetic phosphorimetry method was used to analyze samples containing uranium in the concentration range of 1×10^{-4} $\mu\text{g U/g}$ of solution to 5.0 $\mu\text{g U/g}$ of solution. A Kinetic Phosphorimetry Analyzer KPA-11 (KPA) from Chemchek Instruments Inc. (Richland, WA) was employed for these analyses.¹² Uranium samples for analysis were prepared by dissolving (and diluting if necessary) the sample in a known amount by weight of 1 M HNO₃. The samples were then analyzed according to the procedure specified in the KPA manual.¹³ Two calibration curves were utilized to cover the range of 1×10^{-4} $\mu\text{g U/g}$ solution to 5.0 $\mu\text{g U/g}$ solution. The low curve ranged from 1×10^{-4} $\mu\text{g U/g}$ solution to 1×10^{-2} $\mu\text{g U/g}$ solution and the high curve ranged from 1×10^{-2} $\mu\text{g U/g}$ of solution to 5.0 $\mu\text{g U/g}$ of solution. The accuracy and precision of uranium knowns analyzed using this method was found to be $\pm 5\%$. This accuracy for the method was sufficient to determine elution curves and to estimate uranium losses at various stages in the ion exchange process.

The NBL Titrimetric Method¹⁴ was used to evaluate the quantitative recovery of uranium. The accuracy and precision of this method estimated from the results of 18 uranium control samples analyzed with the test samples was found to be +0.013% with a standard deviation of $\pm 0.047\%$.

X-ray fluorescence spectrometry was used for impurity analyses.¹⁵ Standards for the calibration of the x-ray instrument were prepared from the 25 different atomic spectrometry standard solutions. The 25 elements were combined into three stock solutions (containing 10, 10 and 5 elements, respectively) with a final concentration of approximately 100 μg of each element/g solution. Each stock solution was then diluted successively by a 10:1 weight ratio to yield solutions containing metal ion concentrations of 10 μg element/g solution and 1 μg element/g solution. The separation into 3 different stock solutions was

necessitated by a desire to group elements with similar matrices together to ensure solution stability. The standards were then prepared for x-ray analysis by transferring known quantities of the mixtures onto filter paper and evaporating prior to analysis. These standards, ranging from 0.1 to 80 μg , were used to determine the detection limits of the instrument, and to prepare calibration curves. Known samples prepared and tested against the calibration standards resulted in errors as high as 25% for some elements. The results were used to provide a semi-quantitative picture of the presence of impurities at various stages of the ion exchange process. (In the future, ICP-AES will be used to determine the quantities of impurities separated by the ion exchange in the statistical testing of the method, and the results from that study will be reported separately).

RESULTS AND DISCUSSION

In the initial planning for this project, the following 27 elements were identified as commonly found impurities that could potentially interfere with uranium analysis if present in the sample: Ba, V, Ni, Na, Al, Fe, Ag, Bi, Cu, Co, Mg, Gd, Pr, Cd, Mn, Zn, Ca, Sb, Ti, Cr, Pb, P, Zr, Si, Ta, Sn, and Mo. Procedures capable of separating a majority of these elements from uranium needed to be identified.

The U/TEVA resin was tested because it is very selective for uranium, even in a complex sample matrix. The scheme involves loading and washing the column with 2.5 *M* nitric acid and eluting with 0.01 *M* nitric acid. Recoveries in the 95% range for μg sample sizes have been reported in the literature,⁷ with complete separation from some 30 elemental impurities. The elements separated were Li, Na, Mg, Al, K, Ca, Cr, Mn, Fe, Co, Ni, Cu, Zn, Rb, Sr, Y, Zr, Ru, Rh, Ag, Cd, Cs, Ba, La, Ce, Pr, Nd, Sm, Eu, and Pb. Some variation to this basic procedure and/or sequential treatments can separate V, W, Mo, Bi, Sb, Th(IV) and other actinides.

The investigation into anion exchange resins was prompted by the success of several reported separation methods using them. Based on reported distribution coefficients of uranium and other elements in different media on strongly basic anion resin, the hydrochloric acid media appeared to have the best chance of meeting the requirements of good separation and high uranium recovery. There are several studies related to the adsorption/desorption of uranium and other elements at different HCl molarities.¹⁶

The choice of 8 *M* HCl to adsorb uranium and 1 *M* HCl to desorb it was primarily based on achieving maximum separation from most of the impurities and desorbing all of the uranium. The lanthanides, alkali and alkaline earth metals, Al, Ni, P and Si are not retained on the resin at any acid concentration. Pb(II) is not adsorbed above 6 *M* HCl, while Mn(II) is only adsorbed at or above 9 *M* HCl; therefore these elements are not retained in an 8 *M* HCl wash. While Zn, Cd, Sn(II), Sn(IV), Ta, Bi and Ag are retained in 8 *M* HCl, they are also largely retained in 1 *M* HCl and therefore should remain on the column during uranium elution.

The ions of Ti, Zr, Hf, Mo(VI), Cu(II), and V(V) are potential separation problems. Since they have low distribution coefficients at 8 M HCl acid concentration, it is difficult to predict their actual behavior and the ease of separation of these ions from uranium. This scheme will not separate Fe(III), Mo(VI), Sb(V), or Co(II) from uranium as they are retained under the conditions specified and will be eluted with uranium. As an added benefit, this scheme should readily separate (though it was beyond the scope of this report) a number of actinides from uranium, including Th(IV), Ac(III), Pu(III), and Np(III).

The primary objective of this study was the quantitative recovery of uranium. Thus, it was important to investigate the effects of varying the conditions so as to optimize uranium recovery. Sorption and elution profiles for both the U/TEVA and Bio-Rad anion exchange resins were determined for pure uranium samples. For the Bio-Rad resins only, the methods were optimized and uranium recovery determined with and without the presence of varying levels of impurities.

Sorption and Elution Profiles of Pure Uranium Solutions with U/TEVA

The retention and elution behavior of uranium on U/TEVA resin under four different loading acid concentrations (acid concentrations of 2.5 M HNO₃, 2 M HNO₃, 1.5 M HNO₃ and 1 M HNO₃) was tested. Elution was performed with 0.01 M HNO₃. Samples consisting of 12 mg aliquants of pure uranium were dissolved in 15 mL of acid and the solutions were transferred to preconditioned columns. An additional 10 mL of acid was used to rinse the beaker and was transferred to the column. After passing through the column, the 25 mL loading solution was collected and evaporated. The columns were washed successively with 20 mL and 10 mL aliquots of loading acid. These washings were collected separately. Uranium was then eluted with varying (2 mL to 10 mL) portions of 0.01 M HNO₃ and the eluates were collected to establish elution curves. All the collected fractions were evaporated, fumed twice with concentrated nitric acid, and assayed for uranium using the KPA technique.

The results for sorption and elution of uranium at the various acid concentrations are presented in Table I. Uranium was fully sorbed and strongly retained on the resin in 2.5 M HNO₃. Uranium losses from the sample loading (25 mL volume) and column washings (30 mL volume) totaled less than 1 µg. As the acid concentration of the loading solution was lowered, there was a trend of increased loss of uranium during sample loading and column washes. In 2 M HNO₃, the loss increased to 0.4%, which was unacceptable for our goal of achieving greater than 99.9% recovery. In 1.5 M HNO₃, the loss increased significantly to 6%. In 1 M HNO₃, over 26% of the uranium was lost in the load and washes (55 mL). These results clearly indicated that acid concentrations of 2.0 M or less were unfavorable and were not further evaluated.

Table I Elution Pattern of Uranium on U/TEVA Column					
Volume (mL)	Effluent Fraction	Nitric acid concentration used for Load & Washes			
		1.0 M	1.5 M	2.0 M	2.5 M
		(ug U)	(ug U)	(ug U)	(ug U)
25	Load	8.9	2.8	0.0	0.2
20	Wash 1	1103.1	204.3	7.7	0.1
10	Wash 2	1846.8	474.0	37.8	0.1
02	Eluate 1	3222.8	2001.0	461.5	0.2
03	Eluate 2	5021.5	8773.4	11443.0	3150.6
05	Eluate 3	11.8	13.2	28.8	7084.4
05	Eluate 4	3.3	4.2	5.7	102.2
05	Eluate 5	1.7	2.0	2.8	37.2
05	Eluate 6	0.1	1.6	1.5	21.2
05	Eluate 7	1.0	1.2	1.1	16.1
10	Eluate 8	1.3	1.5	1.7	8.8
Elution was performed with 0.01 M Nitric Acid					

The data in Table I also indicate that a substantial portion of the uranium was still eluting even after addition of 30 mL of eluant. In the 10 mL fraction following the 30 mL eluant addition, all four acid concentrations of the loading solution showed greater than 1 μg eluting. The elution losses versus the acid concentration follow the opposite trend of load and wash losses. As the acid concentration of the loading solution was lowered, tailing during elution decreased. Apparently the higher acid concentration resulted in better retention of uranium during the loading and washing stages, but the uranium was retained to a greater extent during the elution stage and this resulted in significant tailing of the uranium. Based on the combined losses during loading, washing and elution, the maximum achievable recovery of uranium was 98% for the 2.5 M HNO_3 experiment. Further experimentation therefore employed 2.5 M HNO_3 for loading and washing. Because of the significant tailing in the elution stage, the elution characteristics of uranium on U/TEVA were examined using an eluant of 0.1 M ammonium bioxalate buffer.¹⁷ The bioxalate acts to complex uranium and assist in elution, therefore providing greater recoveries. All other experimental conditions were the same as noted above, with the exception that sample loading and washing utilized a 2.5 M HNO_3 solution only. Results are displayed in Table II.

Table II
Sorption and Elution Profiles of U/TEVA resin with 0.01 M Ammonium Bioxalate Eluant

Elution Curve			Uranium Recovery & Elution Data		
Volume (mL)	Sample	Uranium (ug)	Volume (mL)	Sample	Uranium (ug)
25	Load	0.02	25	Load	0.03
30	Wash	0.02	30	Wash	0.06
10	Eluate 1	10685.60	30	DG Analysis	12006.21
10	Eluate 2	14.76	10	Elute 1	2.23
10	Eluate 3	2.52	10	Elute 2	1.15
10	Eluate 4	1.23	10	Elute 3	0.81
10	Eluate 5	0.62	20	Elute 4	0.74
10	Eluate 6	0.92	20	Elute 5	0.83
20	Eluate 7	0.61			
20	Eluate 8	0.82	Average % Recovery:		99.85

-Load and Wash in 2.5 M Nitric
 -Elution with 0.01 M Ammonium Bioxalate
 -All data are averages of 3 trials

While uranium tailing decreased and elution increased with the use of the bioxalate, KPA analysis of the elutions revealed that a significant amount of uranium elution continued well after 40 mL of solution had passed. After the first 40 mL of elution buffer was collected, another 60 mL was added. In the three samples examined, the amount of uranium eluting in this 60 mL was 2.9 μg , 2.6 μg and 3.4 μg , which represented an average loss of 0.025%. The use of ammonium bioxalate as an eluant had another disadvantage. After evaporation of the eluate, large amounts of ammonium bioxalate deposit remained. This salt was somewhat difficult to remove. Heating in concentrated nitric acid reduced the amount of deposit somewhat, and the addition of a few drops of hydrogen peroxide accelerated the decomposition. However, numerous fumings with the nitric-peroxide mixture were necessary to completely remove the salt, and the danger of sample loss due to spattering was significant.

Three samples of pure uranium were recovered by this procedure using 50 mL of the eluant, and were analyzed for uranium recovery using the NBL Titrimetric Method. The mean uranium recovery was 99.85%. Because of the greater recoveries obtained on the initial testing of Bio-Rad ion exchange resins (described below), it was decided to proceed with the Bio-Rad resins only. However, the excellent separation ability of the U/TEVA resin combined with improved recovery with a complexant solution indicates that further experimentation with the U/TEVA resin may be justified in the future.

Sorption and Elution Profiles of Pure Uranium Solutions with Bio-Rad Resins

The performances of four Bio-Rad resins, AG 1-X2, 1-X4, 1-X8 and MP-1, were tested with pure uranium. This allowed an initial comparison of the capabilities of the resins prior to method optimization. Experiments were done in triplicate and repeated to confirm the results. Samples containing 12 mg of uranium were dissolved in 15 mL of 8 M HCl, and the solutions were transferred to preconditioned columns. An additional 10 mL of acid was used to rinse the beaker, and the wash solutions were transferred to the column. The 25 mL of solution used to load the sample was collected. The columns were successively washed three times (20 mL, 10 mL, and 10 mL aliquots) with 8 M HCl and the wash fractions were collected separately. Uranium was then eluted with varying portions (initially with 2 mL and 3 mL aliquots, then five 5 mL aliquots, and finally two 5 mL aliquots were combined to give a 10 mL portion and thus a total of 40 mL of eluant) of 1 M HCl. These eluate fractions were individually collected. All the collected load, wash and eluate fractions were evaporated, fumed twice in concentrated nitric acid, and assayed for uranium using the KPA technique.

Table III gives the sorption and elution behavior of the four resins. When the sample was transferred onto the resin, a pale yellow-colored region was observed at the top of the resin bed. With the AG 1-X2 resin, the KPA data indicated that uranium was lost in both the load and the wash fractions. This was visually observable as the broad, pale yellow band became diffuse and moved down the column. An average of over 30% of the uranium was lost in the load and wash solutions for this resin. The other three resins performed significantly better, with cumulative losses during the 25 mL load and 40 mL wash being less than 5 μ g U. While the majority of the uranium was desorbed in the first 20 mL of eluant, there was some tailing evident in later eluate fractions. Based on these results, the AG 1-X2 resin was eliminated from the study.

Recovery of Pure Uranium with Bio-Rad Resins

Testing was performed on the recovery of uranium on Bio-Rad AG 1-X4, 1-X8 and MP-1 resins using the NBL Titrimetric Method to quantify the recovery of the uranium in the eluate. The method described above, using 8 M HCl with minor modifications (described below), was tested on all four resins. An additional wash step with 5 M HCl was also introduced for testing on each resin. The purpose of this step was to remove several metal ions (Zr, Hf, Cu, Co for example) that have low, yet measurable, distribution ratios in 8 M HCl. It was not certain whether separation of these impurities from uranium would be possible using only 8 M HCl washes. Using 5 M HCl washes, these elements should be readily removed from the resin. However, there was also the possibility of losing uranium, since Uranium has a distribution ratio in 5 M HCl of about 300 compared to approximately 800 in 8 M HCl.

Table III
Adsorption and Elution Behavior of Uranium on Bio-Rad Resins

Volume (mL)	Effluent Fraction	Micrograms Uranium in Fraction			
		Bio-Rad Resin			
		1-X2	1-X4	1-X8	MP-1
25	Load	21.0	0.3	0.5	0.1
20	Wash 1	740.1	0.6	0.5	0.1
10	Wash 2	1131.6	1.4	0.3	0.7
10	Wash 3	1754.3	2.5	0.1	4.0
2	Eluate 1	7761.9	8.0	644.9	0.3
3	Eluate 2	81.6	9719.8	7646.7	9918.4
5	Eluate 3	7.2	1484.1	2152.7	890.9
5	Eluate 4	0.4	13.8	145.4	110.2
5	Eluate 5	0.1	2.5	17.7	14.7
5	Eluate 6	0.0	0.5	3.5	3.0
5	Eluate 7	0.0	0.1	0.8	1.5
10	Eluate 8	0.0	1.8	0.3	0.2

-Load and Washes done in 8 M HCl

-Elutions in 1 M HCl

-Data are averages of 3 runs with the exception of 1-X2 resin (2 runs)

All samples were loaded onto columns in the same manner used for the elution profiles. For one set of samples, the columns were washed by three 10 mL additions of 8 M HCl (30 mL total volume). For the second set of samples, the columns were washed with three 5 mL additions (15 mL total volume) of 8 M HCl. This was followed by three 5 mL additions (15 mL volume) of 5 M HCl. All fractions were collected, evaporated and fumed for analysis by the KPA method. Uranium was then eluted from all samples with three 10 mL portions (30 mL total volume) of 1 M HCl for analysis by the NBL Titrimetric Method. An additional 10 mL portion of 1 M HCl was added and collected for assay by the KPA method. Sample beakers were also tested for any un-transferred uranium by the KPA method. Triplicate analyses were performed for all resin and washing schemes. Results are presented in Table IV.

Table IV Uranium Recovery with Bio-Rad Resins							
Volume (mL)	Effluent Fraction	Bio-Rad Resin					
		1-X4		1-X8		MP-1	
		8M Only	8M & 5M	8M Only	8M & 5M	8M Only	8M & 5M
		(ug U)	(ug U)	(ug U)	(ug U)	(ug U)	(ug U)
25	Load	0.2	0.4	0.2	0.1	0.1	0.2
30 or 15	Wash (8 M)	18.0	0.6	0.4	0.2	0.3	3.1
15	Wash (5 M)	-	108.8	-	2.2	-	15.5
30	Eluate (DG)	12103.6	11502.8	12121.1	12014.7	12006.3	11970.0
10	Eluate 2	6.1	2.8	16.7	11.4	4.2	19.9
% Recovery:		99.75	98.81	99.79	99.82	99.85	99.63
-Recovery determined by the NBL Titrimetric Method -Elution performed with 1 M HCl -Combined 8M and 5M HCl washes performed with 15 mL each							

For all three resins, using 8 M HCl washes yielded an average uranium recovery of 99.8%. The AG 1-X4 resin, however, exhibited larger uranium losses during the wash steps than the other two resins. In order to minimize these losses, a reduction in load and wash volumes would probably be necessary. However, that would compromise loading efficiency and increase the probability of unadsorbed and slightly adsorbed impurities remaining on the column prior to uranium elution. The data also shows a significant loss of uranium, much higher than seen for the other resins, upon washing the AG 1-X4 resin with 5 M HCl. This was not unexpected due to the lower distribution coefficients for uranium on lower cross-linked resins. For these reasons, the AG 1-X4 resin was eliminated from further experimentation.

There were no significant uranium losses with the AG 1-X8 and AG MP-1 resins during the loading and washings, and recoveries with AG 1-X8 and AG MP-1 averaged 99.8% in the 30 mL of eluant collected, with some uranium found in the next 10 mL of the eluate collected. Though the recoveries in the AG 1-X8 and AG MP-1 resins were approaching the goal, the presence of uranium beyond the 30 mL of eluate was unexpected. Our initial elution behavior studies with the resins involving only the 8 M HCl washes indicated very little uranium tailing should be expected following 30 mL of elution. However, significant amounts of uranium were now detected by KPA in the final 10 mL of eluate.

One significant difference in the elution behavior studies and the uranium recovery experiments involved the addition of the eluant. During elution behavior studies, eluant was added in 5 mL portions or less. For the uranium recovery test, eluant was added in 10 mL portions. While the total amount of eluant used was the same, it was speculated that adding more aliquots of smaller volumes improved the overall exchange kinetics, thus minimizing tailing of the uranium.

To optimize uranium recovery, uranium was eluted using eight 5 mL portions of 1 M HCl for a total of 40 mL. These conditions were tested in duplicate on the AG 1-X8 and AG MP-1 resins using both an 8 M and a combination 8 M/5 M HCl wash as before. All other experimental details remained unchanged. The results are given in Table V. The NBL Titrimetric Method results indicated that uranium recovery for both the AG 1-X8 and AG MP-1 resins averaged 100.0%. KPA analysis corroborated these results, indicating no significant losses in the load or wash, or in the 10 mL of eluant following the first 40 mL of eluate. With this experiment, the ion exchange procedure to be used for further testing with impurities was finalized.

Table V Optimized Uranium Recovery with Bio-Rad 1-X8 and MP-1 Resins					
Volume (mL)	Effluent Fraction	Bio-Rad Resin			
		1-X8		MP-1	
		8M Only	8M & 5M	8M Only	8M & 5M
		(ug U)	(ug U)	(ug U)	(ug U)
25	Load	0.1	0.2	0.0	0.0
30 or 15	Wash (8 M)	0.1	0.1	0.1	0.0
15	Wash (5 M)	-	2.1	-	3.5
40	Eluate (DG)	11935.0	11937.7	11916.9	12036.4
10	Eluate	0.3	0.1	0.4	0.0
% Recovery:		100.01	100.01	100.06	99.87
-Recovery determined by the NBL Titrimetric Method -Elution performed with 1 M HCl -Combined 8M and 5M washes performed with 15 mL each					

Ion Exchange of the Matrix Blank

The optimized method described above was used to study the adsorption and elution behavior of elements in the matrix blank (no uranium present). Three matrix blanks containing 50 µg of each of the 25 elements and one reagent blank (no uranium or impurities) were ion exchanged following the method described

above. The Ag and Mo impurities were removed from the original list of 27 elements due to the incompatibility of their standard solutions with the other impurities. A reagent blank was run with the matrix blanks to take into account any contributions from the reagents and/or any contamination. From each procedure, five portions were collected: the load, the 8 M HCl wash, the 5 M HCl wash, the 1 M HCl eluate and a 1 M HCl wash of the empty sample beaker. The first four portions were analyzed by x-ray fluorescence spectrometry. The 1 M HCl wash of each empty sample beaker was analyzed for trace uranium by the KPA method.

Results obtained on the ion exchange of matrix blanks are given in Table VI. The results of only 22 elements are given in this table. Sb, Sn and Cd are not included for reasons explained below. All of the theoretically predicted non-complexing metal ions (alkali, alkaline earth, Ni, Pb, Mn) were eluted in the load solutions. Though the results were difficult to interpret, it appears that aluminum was also removed in the loading step as expected. It was uncertain whether elements with low but appreciable distribution coefficients, such as Cu(II), Ti, Zr, Co and V(V), would be removed from the resin with 8 M HCl. However, from the 8 M HCl load results, it appears that Ti, Zr, V(V), and most of the Cu were removed in the 8 M HCl washes, and only Co eluted in the 5 M HCl wash step.

The 1 M HCl eluate contained predominately Fe, as expected. While Sb was also expected to elute with uranium, it was difficult to measure Sb (as well as Sn and Cd) with any certainty due to the proximity of their primary lines with the Re lines resulting from the Re filament of the instrument's x-ray tube. The elements Bi, Zn, Cd, Ta and Sn that were expected to adhere strongly to the resin at 1 M HCl did not appear in any effluent and therefore are assumed to remain adsorbed to the resin bed. It became clear with these results that 5 M HCl washes were not necessary since most of the metal impurities, with the exception of Co, can be separated from uranium with 8 M HCl sample adsorption solution and 1 M HCl elution solution.

Ion Exchange of Spiked Uranium Samples

The optimized method was tested on 12 mg uranium samples spiked with approximately 50 μg of each of the 25 elements in order to determine how the method performed on impure samples. Both the AG 1-X8 and MP-1 resins were tested in duplicate with 8 M HCl washes only and sequential 8 M HCl and 5 M HCl washes. The load and washes were analyzed by x-ray fluorescence for the impurities to confirm the behavior of these elements in the presence of uranium, and the eluate was tested for uranium recovery by the NBL Titrimetric Method. A 1 M HCl wash of each empty sample beaker was analyzed for trace uranium by the KPA method.

Table VI
Sorption and Elution Pattern of 22 Impurities on Bio-Rad MP-1 Resin

Element	8 M HCl LOAD (25 mL)		8 M HCl WASH (15 mL)		5 M HCl WASH (15 mL)		1M HCl ELUATE SOLUTION (40 mL)	
	Blank	Sample	Blank	Sample	Blank	Sample	Blank	Sample
	ug Impurity	ug Impurity	ug Impurity	ug Impurity	ug Impurity	ug Impurity	ug Impurity	ug Impurity
Ba	0	41.9	0	0	0	0	0	0
Na	3.6	18.4	2.4	4.7	16.1	3.1	3.9	4.5
Bi	0.7	2	3.2	3.4	3.5	3	3.6	4
Ni	0.7	32	0.5	1.7	1	1	0.5	0.5
V	1.1	34.7	1.1	1.4	1.1	1.1	1.1	1.1
Fe	3.1	1.7	1.8	5.3	1.1	3.1	3	49
Al	0	9.3	0	0	0	0	0	0
Cu	0.8	6.1	1.4	23.6	1.3	8.3	1.2	1.6
Co	0.6	0.8	0.9	5.4	1	73.5	1	1
Mg	0.9	9.9	1	1.1	1	0.9	0.9	0.9
Gd	0	37.5	0	2.2	0	0	0	0
Pr	0	37.1	0	0	0	0	0	0
Zn	0	0	0	0.5	0.1	0.7	0	0
Ca	2.9	37.9	3.7	3.6	3.9	1.3	0.8	0
Ti	0	-	0	0.9	0	0	0	0
Cr	0	36.2	0	1.4	0	0	0	0
Pb	1.7	36.3	4.5	6.6	4.5	4.8	4.3	4.5
Mn	0	37.7	0	2.4	0	0	0	0
P	0.1	31.5	0.1	1.2	0.2	0.2	0.1	0.1
Zr	0	54.8	2.4	3.3	1.1	1.1	0.9	1.6
Si	16.1	92.1	3.9	15.4	5.5	5.4	24.9	4.4
Ta	0	1.3	0	2.2	2.4	1	2.4	0

-Results determined by x-ray fluorescence. Measurement uncertainties averaged +/- 3.0 ug, though several uncertainties were as high as +/- 12 ug.

The results from the analysis of the impurities showed no changes from those obtained with the matrix blank. The elements were eluted or retained, as before. However, uranium recoveries were found to vary from 99.5% to 99.9%, with an average recovery of 99.8%. To investigate this loss of uranium, the load and wash solutions were checked qualitatively for uranium by x-ray fluorescence spectrometry. In some cases uranium was detected in load fractions and/or the corresponding wash fractions. Two possible causes for the loss were identified: an overload of the capacity of the column, or the presence of an interfering matrix element or elements. Since the total amount of sample loaded, uranium and impurities, does not exceed 10% of the column capacity, it appeared that one or more of the impurities present was causing the breakthrough of uranium.

The role of the impurities on the uranium loss was investigated. The elements that are known to form complexes or interfere with uranium were identified. These elements included P, Si, and F. The F was not specifically added to the uranium spike solutions, but some of the metal impurities were stabilized in solution as fluoride complexes. In addition, Sb, Fe, and Co were tested for interference with uranium recovery. Antimony was selected because it was added as a tartrate solution, and tartrate may complex uranium. Iron and Co (among others) are adsorbed onto the resin along with uranium and thus reduce the number of resin sites available for uranium adsorption.

The method was tested for recovery on samples containing 12 mg uranium mixed with selected impurities, each at the 50 μg level (unless impurity level is stated otherwise). All tests were done in duplicate. Description of samples ion exchanged and tested for uranium recovery are as follows: uranium and added P; uranium and added Si; uranium and added Na; uranium and added Sb (with tartrate); uranium and all impurities except P; uranium and all impurities except P and Si; uranium and all impurities except P, Si, Ti, Sn, Ta, Zr (no F present); uranium and Sb (with tartrate); uranium and all impurities but Sb (and tartrate); uranium and all impurities but Co and Sb (and tartrate); uranium and all added impurities at the 15 μg level each; and finally uranium and all added impurities at the 30 μg level each. The results of these experiments are found in Table VII.

Uranium recoveries from samples with P and Si individually show 100.0% recovery. The results of samples with uranium and all impurities except for P, Si, Ti, Sn, Ta, Zr (no F present) was 99.8%. Therefore P, Si and F were ruled out as the cause of uranium losses.

The results of uranium with impurities minus Sb (with tartrate), and uranium with impurities minus Sb (and tartrate), Fe and Co resulted in 100.0% recovery. Based on this result, Sb (with tartrate) was initially suspected to cause uranium losses. However, uranium recovery with uranium plus Sb (with tartrate) was 100.0%. Samples of uranium with all impurities (including Sb with tartrate) at the 15 μg and 30 μg levels each resulted in 100.0% recoveries. While it was not shown that Sb (with tartrate) was the cause of the uranium loss, it was removed from the study due to the uncertainty of the effect of tartrate on the ion exchange

behavior of both uranium and antimony. In order to test the method's ability to handle greater impurity levels, spiked 12 mg uranium samples were tested in triplicate at the 50 μg , 100 μg and 200 μg level each with 24 impurities (no Sb with tartrate). The results, shown in Table VII, indicate 99.98% uranium recoveries.

Table VII
Uranium Recovery of Spiked Samples

Sample Description	Resin	% U Recovery
U + Imp Group A	1-X8	99.85
U + Imp Group A; 8M & 5M HCl Wash	1-X8	99.70
U + Imp Group A	MP-1	99.88
U + Imp Group A; 8M & 5M HCl Wash	MP-1	99.82
U + P	1-X8	100.02
U + Si	1-X8	100.09
U + Imp Group A - P	1-X8	99.88
U + Imp Group A - P - Si	1-X8	99.85
U + Imp Group B - P - Si - F	1-X8	99.82
U + Sb	1-X8	99.98
U + Imp Group A - Sb	1-X8	99.99
U + Imp Group A - Sb - Co - Fe	1-X8	100.02
U + Imp Group A (15 ug each)	1-X8	100.01
U + Imp Group A (30 ug each)	1-X8	100.02
Optimized Method (no Sb)		
U + Imp Group A - Sb (50 ug each)	1-X8	99.98
U + Imp Group A - Sb (100 ug each)	1-X8	99.98
U + Imp Group A - Sb (200 ug each)	1-X8	99.97
<ul style="list-style-type: none"> - "U + Imp" is uranium with all impurities added - Recoveries determined by the NBL Modified Titrimetric Method - Imp Group A is Ba, Na, Bi, Ni, V, Fe, Al, Cu, Co, Mg, Gd, Pr, Zn, Ca, Ti, Cr, Pb, Mn, P, Zr, Si, Ta, Sn, Sb, Cd. - Imp Group B is Imp Group A minus Ti, Sn, Ta and Zr (no F present) 		

Uranium Breakthrough and Recovery

The Bio-Rad AG 1-X8 and AG MP-1 resins were tested for uranium breakthrough and uranium recovery on 5 mg, 10 mg, 15 mg, 20 mg and 25 mg size uranium samples. The method involving the 5 *M* HCl washes was eliminated from further testing since it added little additional selectivity.

For uranium breakthrough tests, samples were transferred with 25 mL of 8 *M* HCl and successively washed with 20 mL of 8 *M* HCl until at least 160 mL of the acid was added to the column. The fractions were collected, dried, fumed twice with concentrated nitric acid and then the uranium analyzed by KPA.

For uranium recovery tests, samples were transferred with 25 mL of 8 *M* HCl, washed with four 5 mL portions of 8 *M* HCl, and eluted with eight 5 mL portions of 1 *M* HCl. The three effluents were evaporated and fumed twice with concentrated nitric acid. The load and wash samples were analyzed by KPA to check for uranium losses and the eluate was analyzed by the NBL Titrimetric Method for uranium recoveries.

The results, shown in Table VIII, for similar sample sizes indicate that there is more uranium breakthrough with the AG MP-1 resin than the AG 1-X8 resin for similar size samples. The total amount of all uranium in the 160 mL washes was found to be as high as 15 μg on the 5 mg and 10 mg samples. In comparison the actual losses with higher size samples did not exceed 5 μg uranium. How the difference in sample sizes affects uranium breakthrough is not clear; however, the differences in AG 1-X8 and AG MP-1 behavior can be explained on the basis of resin fines present in AG MP-1, which may be responsible for uranium breakthrough. The overall results indicate that the losses ($< 0.02\%$) are minimal for the amounts of acid (< 50 mL) used for the ion exchange of actual samples.

Statistical Testing

In order to obtain a more valid determination of the total uranium recovery from the ion exchange method described in this report and to validate this method, a statistical plan was designed. Two resins, Bio-Rad AG 1-X8 and AG MP-1, were tested using the 8 *M* HCl wash only. The U/TEVA method was not tested at this time since the initial recovery data did not achieve our goal of greater than 99.9% uranium recovery.

Samples were loaded with 25 mL of 8 *M* HCl, washed with 20 mL of 8 *M* HCl, and eluted with 40 mL of 1 *M* HCl. Samples were aliquanted into glass and Teflon beakers. The samples ion exchanged were: reagent blanks, 50 μg each of 24 impurity elements as matrix blanks, 12 mg uranium, and 12 mg uranium spiked with 50 μg each of 24 impurity elements.

The 24 impurity elements used in this statistical plan were Ba, Na, Bi, Ni, V, Fe, Al, Cu, Co, Mg, Gd, Pr, Zn, Ca, Ti, Cr, Pb, Mn, P, Zr, Si, Ta, Sn, and Cd.

Table VIII
Uranium Breakthrough Study on Bio-Rad 1-X8 and MP-1 Resins

MP-1 Resin					
Wash Step	Uranium Sample Size				
	5 mg	10 mg	15 mg	20 mg	25 mg
1	0.01	0.04	0.00	0.00	0.00
2	0.71	1.32	0.10	0.05	0.06
3	2.22	3.79	0.66	0.09	0.31
4	1.92	3.38	1.14	0.99	0.53
5	1.07	2.24	0.92	1.41	0.83
6	0.60	1.31	0.92	0.95	0.94
7	0.25	0.55	0.59	0.55	0.30
8	0.38	0.78	0.58	0.30	0.02
1-X8 Resin					
Wash Step	Uranium Sample Size				
	5 mg	10 mg	15 mg	20 mg	25 mg
1	0.40	1.11	0.09	0.10	0.65
2	0.15	0.88	0.16	0.16	0.28
3	0.27	0.44	0.14	0.24	0.29
4	0.17	0.30	0.17	0.18	0.31
5	0.21	0.34	0.07	0.10	0.24
6	0.09	0.06	0.17	0.04	0.17
7	0.12	0.04	0.15	0.04	0.10
8	0.10	0.14	0.04	0.20	0.07
-Results are ug U found in each wash step -Each wash was performed with 20 mL of 8 M HCl -Analysis of eluate performed by KPA analysis					

The statistical plan and data analysis were provided by an NBL statistician. The scheme was designed to provide information on the recovery of uranium with and without impurities, the day-to-day variation in ion-exchange, the analyst-to-analyst variation, the recovery of non-ion exchanged versus ion exchanged samples, and the performance of glass versus Teflon beakers. Reagent blanks and matrix blanks were included in the statistical plan, but are not included in this report. All samples were aliquanted by one

analyst, and all NBL Titrimetric Method titrations were performed by one analyst. The ion exchange procedure was performed by two analysts, each processing 10 samples each day for 4 days.

During each day of titration, three "blind" control samples were analyzed by the analyst performing the NBL Titrimetric Method. Of the 9 "blind" controls run during the titrations of AG 1-X8 resin separated samples, the mean was 100.031% of the expected value with standard deviation of $\pm 0.055\%$. For the 9 "blind" controls run during the titrations of AG MP-1 resin separated samples, the mean was 99.996% of the expected value with standard deviation of $\pm 0.031\%$.

Table IX Summary Table of Resin Statistics for Uranium Recovery						
	BIO-RAD AG 1-X8			BIO-RAD MP-1		
	n	Mean % RD	Std Dev	n	Mean % RD	Std Dev
Overall:	30	-0.073	0.085	27	-0.045	0.075
ANALYST AND DAY						
Analyst 1:						
Day 1	7	0.004	0.034	7	-0.009	0.075
Day 2	8	-0.045	0.068	7	-0.053	0.073
Total:	15	-0.022	0.059	14	-0.031	0.075
Analyst 2:						
Day 1	8	-0.089	0.044	7	-0.030	0.049
Day 2	7	-0.162	0.095	6	-0.095	0.089
Total:	15	-0.123	0.079	13	-0.060	0.075
SAMPLE TYPE						
Non-Ion exchanged:	9	0.012	0.043	9	-0.034	0.052
Ion exchanged:						
Pure	16	-0.092	0.101	13	-0.028	0.068
Impurities added	14	-0.051	0.060	14	-0.060	0.081
D&G TITRATION DAY						
Day 1	8	-0.118	0.107	6	-0.083	0.083
Day 2	10	-0.082	0.087	10	-0.042	0.077
Day 3	12	-0.043	0.058	11	-0.027	0.068
BEAKER TYPE						
Glass	15	-0.072	0.086	14	-0.045	0.059
Teflon	15	-0.074	0.088	13	-0.045	0.091
PAIRED DIFFERENCE						
RDno Imp - RDimp	14	-0.041	0.078	13	0.036	0.090

Results from the statistical analysis of the data are summarized in Table IX. Both resins exhibited similar recoveries and standard deviations for samples spiked with impurities. The mean recovery for samples separated by the AG 1-X8 resin was 99.95%, and for samples separated by the AG MP-1 resin it was 99.94%. Standard deviations for the two were $\pm 0.060\%$ and $\pm 0.081\%$, respectively. There were some differences in the recoveries of the samples containing only pure uranium. The mean recovery of samples separated by the AG 1-X8 resin was 99.91%, and for samples separated by the AG MP-1 resin it was 99.97%. The standard deviations also differed somewhat and were $\pm 0.101\%$ and $\pm 0.068\%$, respectively. The mean recovery of uranium samples (not ion exchanged) assayed with the AG 1-X8 ion exchanged samples was 100.01%, and for uranium samples (not ion exchanged) assayed with the AG MP-1 ion exchanged samples it was 99.97%. The standard deviations for these samples were very similar and were $\pm 0.043\%$ and $\pm 0.052\%$, respectively.

There was no significant day-to-day titration variation observed, and no difference was observed between samples ion exchanged using glass or teflon beakers. However, there were analyst-to-analyst and ion-exchange day-to-day differences observed. Considering the overall statistical outcome it can be stated that the separation method utilizing either AG 1-X8 or AG MP-1 resin will yield in greater than 99.9% recovery with standard deviations of less than $\pm 0.1\%$.

CONCLUSIONS

This study has shown that quantitative recovery ($>99.9\%+$) of uranium in the presence of numerous impurities can be achieved using the anion exchange procedure outlined above. Both Bio-Rad AG 1-X8 and AG MP-1 resins are capable of achieving the desired recovery while providing separation from the impurities. Future work will concentrate on further improving the recovery and precision of the method, along with determining the uranium sample size range applicable to the method. In addition, quantification of the impurities present in the loads and washes, along with any trace impurities present in the uranium eluate fraction, will be performed using an ICP-AES method.

REFERENCES

1. J. Korkisch, "Handbook of Ion Exchange Resins: Their Application to Inorganic Analytical Chemistry," Volume I and II, CRC Press, Boca Raton, FL, 1989.
2. I. M. Kolthoff, and P. J. Elving, "Treatise on Analytical Chemistry," Part II, Vol. 9, New York, Interscience, 1962.
3. C. E. Pietri, B. P. Freeman, and J. R. Weiss, "The Quantitative Ion Exchange Separation of Plutonium from Impurities," NBL-298 (September 1981).

4. W. G. Mitchell, M. I. Spaletto, K. Lewis, M. D. Soriano, and M. M. Smith, "The Effect of Ion Exchange Purification on the Determination of Plutonium at the New Brunswick Laboratory," NBL-323 (July 1990).
5. H. O. Haug, "Gmelin Handbook of Inorganic Chemistry, Uranium Supplement, D3, Anion Exchange," Springer-Verlag, Berlin, 1982.
6. F. W. E. Strelow and C. H. S. W. Weinert, "Selective Separation of Uranium from Other Elements by Cation Exchange Chromatography in Hydrobromic Acid and Hydrochloric Acid-Acetone Mixtures," *Talanta* Vol. 20, (1973) pp. 1127-1137.
7. E. P. Horowitz, M. L. Dietz, R. Chiarizia, H. Diamond, A. M. Essling, and D. Graczyk, "Separation and Preconcentration of Uranium from Acidic Media by Extraction Chromatography," *Analytica Chimica Acta*, 266 (1992) pp. 25-37.
8. P. M. Santoliquido, "Preparation and Provisional Certification of NBL Spectrographic Impurity Standards CRM 123 (1-7) and CRM 124 (1-7)," NBL-310, September 1983.
9. "Purification of Uranium for Mass Spectrometric Analysis using U/TEVA SPEC Columns," NBL-SP-U(I)-4, NBL Procedures Manual.
10. "Preparation of Uranium Standard Solution," NBL-CAL-U(E)-2, NBL Procedures Manual.
11. "Guide to Ion Exchange," BIO-RAD, Catalog Number 140-9997.
12. R. Brina and A. G. Miller, "Direct Detection of Trace Levels of Uranium by Laser-Induced Kinetic Phosphorimetry," *Anal. Chem.*, Vol. 64, No. 13, 1992, pp. 1413-1418.
13. KPA-11 Operation and Service Manual, Kinetic Phosphorescence Analyzer, Chemchek Instruments, Inc., Richland, WA, June, 1992.
14. "Determination of Uranium by Ferrous Reduction in Phosphoric Acid and Titration with Dichromate (NBL Titrimetric Method)," NBL-SA-U(E)-1, NBL Procedures Manual.
15. Instruction Manual, Sequential X-ray Spectrometer Systems, Rigaku Industrial Corporation, Osaka, Japan.
16. T. N. van der Walt, F. W. E. Strelow and R. Verheij, "The Influence of Crosslinkge on the Distribution Coefficients and Anion Exchange Behaviour of Some Elements in Hydrochloric Acid," *Solvent Extraction and Ion Exchange*, 3(5), (1985) pp. 723-740.
17. EIChroM Industries Inc., SpecNews Volume 1, Issue 1, June 1992.

APPENDIX A

PROCEDURE FOR THE PURIFICATION OF URANIUM BY ANION EXCHANGE

PURIFICATION OF URANIUM BY ANION EXCHANGE SEPARATION (MINI COLUMN)

I. INTRODUCTION

A. Applicability

The anion exchange purification of uranium removes various impurities in order to permit reliable determinations of uranium content by various techniques, such as coulometry and titrimetric methods. This procedure may also be used to separate impurities for impurity analysis by removing all the uranium from the sample. In order to achieve the stated quantitative recovery, uranium sample sizes of 12 mg and impurities not to exceed 10% of total sample weight are optimal.

B. Summary of Method

Uranium, in sample aliquants previously evaporated, is dissolved in 8 M HCl. The U(VI)-chloro anionic complex formed in 8 M HCl is sorbed on Bio-Rad AG 1-X8 or Bio-Rad AG MP-1 anion exchange resin. Impurities are separated by washing the resin with 8 M HCl. Uranium is desorbed from the resin with 1 M HCl. The eluted uranium samples are then dried down. Concentrated nitric acid, and if needed, sulfuric acid, are added and the sample is fumed to dryness.

C. Method Performance

A study was undertaken by Narayanan, et. al.¹ to develop an ion exchange method for removing major impurities and achieving complete uranium recovery. The study included evaluation of the quantity of uranium lost during the ion exchange process. High-purity uranium samples were analyzed with and without prior ion exchange. Also, uranium samples were ion exchanged with and without added impurities. The mean difference in uranium concentration between ion exchanged and non-ion exchanged aliquants of pure uranium samples was found to be -0.104% for Bio-Rad AG 1-X8 resin and +0.006% for Bio-Rad AG MP-1 resin. The mean difference in uranium concentration between pure uranium and uranium with impurities samples was found to be -0.041% for Bio-Rad AG 1-X8 resin and +0.032% for Bio-Rad AG MP-1 resin.

D. Interferences/Limitations

Each batch of anion exchange resin must be tested prior to use to verify that it will provide the stated quantitative recovery.

The choice of 8 M HCl to adsorb uranium and 1 M HCl to desorb was based on achieving maximum separation from most of the impurities while desorbing all of the uranium. The lanthanides, alkali and alkaline earth metals, Al, Ni, P and Si are not retained on the resin at any acid concentration. Pb(II) is not adsorbed above 6 M HCl, while Mn(II) is only adsorbed at or above 9 M HCl; therefore these elements are not retained on the column after an 8 M HCl wash. While Zn, Cd, Sn(II), Sn(IV), Ta, Bi and Ag are retained in 8 M HCl, they are also largely retained in 1 M HCl and therefore should remain on the column during uranium elution. The ions of Ti, Zr, Hf, Mo(VI), Cu(II), and V(V) have low distribution coefficients at 8 M HCl acid concentration. It is therefore difficult to predict their actual behavior and the ease of separation of these ions from uranium. This scheme will not separate Fe(III), Mo(VI), Sb(V), or Co(II) from uranium, as they are retained under the conditions specified and will be eluted with uranium. As an added benefit, this scheme should readily separate a number of actinides from uranium, including Th(IV), Ac(III), Pu(III), and Np(III).

E. Traceability to Primary Standards

This ion exchange method is a purification performed prior to quantitative determination of elemental concentration, isotopic analysis or impurity analysis; traceability to primary standards is described in the analytical procedure used following purification.

To ensure that the ion-exchange separation is performed quantitatively, "blind" standards are ion exchanged and analyzed with each group of samples on which the uranium concentration is being determined. These "blind" standards are aliquants of solutions prepared from CRM 112-A, Uranium Metal Assay Standard, or other appropriate uranium standards, which contain amounts of uranium unknown to the analyst. Results obtained on the blinds are entered into the NBL computer-based quality control program. Poor ion-exchange recoveries will lead to out-of-control conditions.

F. Safety Procedures

Review Material Safety Data Sheets (MSDS) for corrosive acids prior to performing the procedure.

All handling of uranium solutions and uranium aliquants is performed in approved laboratories. Fuming of uranium materials is done on steam baths or on hot plates in a filtered hood. The requirements for working with radioactive materials, such as uranium, may be found in the NBL Radiological Control (RadCon) Manual.

Nitric, sulfuric, and hydrochloric acids are extremely corrosive; gloves should be used when preparing solutions of these acids. Polyethylene or Teflon bottles used to store 8 M HCl should be regularly inspected for hairline cracks before refilling.

II. INSTRUMENTS, APPARATUS, AND REAGENTS

A. Instruments

Not applicable

B. Apparatus

1. Pipets, as needed. One piece disposable plastic transfer pipets, approximate draw of 3.5 mL.
2. Columns, Bodman Polypropylene Columns, Large Capacity Robotic Tube, Order No. 3300-25; polyethylene 70 μ frits, Order No. 3201-00, from Bodman Industries, P. O. Box 2421, Aston, PA 19014.
3. Rod, glass or plastic.
4. Column Holder, NBL design,² preferably painted to avoid corrosion from HCl.
5. Hot Plate.
6. Bottles, polyethylene or Teflon™, 250 mL and 1 L.
7. Tray, enamel-coated metal.
8. Parafilm,® Grade M.

9. Glass Wool.
10. Beakers, glass or Teflon™, 50 mL and 100 mL.

NOTE: If new glassware is used, clean by heating for 2 hours in approximately 8 *M* nitric acid, then rinse with deionized distilled water and finally allow to air dry. To clean used glassware, soak first in approximately 8 *M* nitric acid, then rinse with deionized distilled water. Rinse the inside of the beakers twice with approximately 8 *M* nitric acid followed by rinsing three times with deionized distilled water.

11. Marker pen.

C. Reagents

Reagents are ACS reagent grade unless otherwise specified. Solutions may be prepared in multiples or fractions as required to minimize waste.

1. Water, deionized distilled.
2. Resin, Bio-Rad AG® 1-X8 or Bio-Rad AG® MP-1, 100-200 mesh (dry) from Bio-Rad.
3. Hydrochloric acid, 8 *M* HCl. Add 660 mL of concentrated HCl to about 300 mL of deionized distilled water; dilute to 1 L. Store the solution in a polyethylene or Teflon bottle.
4. Hydrochloric acid, 5 *M* HCl. Add 413 mL of concentrated HCl to about 550 mL of deionized distilled water; dilute to 1 L. Store the solution in a polyethylene or Teflon bottle.
5. Hydrochloric acid, 1 *M* HCl. Add 82 mL of concentrated HCl to about 900 mL of deionized distilled water; dilute to 1 L. Store the solution in a polyethylene or Teflon bottle.
6. Nitric acid, concentrated.
7. Nitric acid, 8 *M*. Add 500 mL concentrated HNO₃ to about 400 mL of deionized distilled water; dilute to 1 L. Store the solution in a Teflon bottle.

8. Sulfuric acid, 3 *M*. Add 168 mL of concentrated sulfuric acid (Baker Instra-Analyzed® or equivalent) to about 800 mL of deionized distilled water, and dilute to 1 L. Store the solution in a polyethylene or Teflon bottle.

III. PROCEDURE

A. Column Preparation

NOTE: The uranium ion exchange is described for a single aliquant. However, standard practice is to perform ion-exchange separations in groups of ten. Follow the same procedure for each aliquant.

1. Clean the bench top inside the hood. Clean the column holder with a wet paper towel. Place the column holder on a clean sheet of paper towel. Place a clean paper towel on the column holder base.
2. Place the frit firmly inside an empty column by gently pushing it to the end of the column with a glass or plastic rod.
3. Place the ion-exchange column in the column holder. Center a 100-mL beaker under the column. The column tip should not touch the sides of the beaker or the liquid which will be collected in it (about 60 mL).
4. Add about 3 mL of deionized distilled water into the column. If no flow is observed in 30 seconds, apply pressure by covering the column opening with the palm of your hand and pressing downward. This will wet the glass frit and start the flow.
5. Prepare the resin in a clean 250-mL polyethylene bottle by mixing dry resin with enough deionized distilled water to make a slurry.
6. Add enough resin to the column by using a disposable pipet to give a resin height (after the water has stopped dripping) of 4 cm measured from the top of the frit in the bottom of the column to the top of the resin. This equals not less than 2.5 mL of wet resin volume. Retain pipet for use in steps III.A.7 and III.A.9.

7. Using the same disposable pipet from step 6 above, add about 5 mL of deionized distilled water to the column. Gently draw some resin from the column into the disposable pipet. Release the resin from the pipet and allow it to settle back into the column. This will help to eliminate air bubbles in the resin and prepare a uniform resin bed.
 8. Use a marker pen to mark 2-mL and 1-mL solution heights on a new disposable pipet. This pipet will be used for 8 M HCl reagent additions to the column.
 9. Add 10 mL of 8 M HCl to the column using the marked disposable pipet. Using the disposable pipet from step 7, again draw some resin from the column into the pipet and release it to gently loosen the resin bed. Allow the bed to settle. Allow solution to drain to top of resin bed.
- NOTE: When the upper liquid level drains to the top of the resin, the column will stop dripping but the resin bed will not go dry. Therefore it is not necessary to add more solution to the column before it stops flowing.
10. For only the columns filled with the Bio-Rad AG MP-1 resin, position a piece (approximately 2 cm in diameter) of glass wool at the top of the resin bed using the plastic or glass rod.
 11. Add 10 mL of 8 M HCl to the column using the marked disposable pipet. Repeat.

B. Sample Preparation

NOTE: Use the same marked disposable pipet used in steps 8 and 9 for reagent additions to the sample.

1. Remove the Parafilm beaker cover from the U sample beaker and invert the Parafilm over the sample beaker.
2. Rinse the surface of the Parafilm into the beaker containing the sample with 2 mL of 8 M HCl. Wash down the beaker walls with an additional 12 to 14 mL of 8 M HCl. Swirl gently until the sample dissolves.

NOTE: Gentle warming of the solution on a hot plate for 5-10 minutes may hasten dissolution of difficult-to-dissolve samples. Cool samples to room temperature before proceeding.

C. Sample Loading

NOTE: If impurities are to be collected for analysis, place a clean 50-mL beaker under the column before beginning sample loading. However, even if not collecting impurities it is a good practice to collect and examine the sample load and rinse solutions for uranium losses.

1. Using the marked pipet, add 2 mL of 8 *M* HCl to the reservoir and let it drain to the top of bed.
2. Wet the sample beaker spout by adding a drop of 8 *M* HCl with the disposable pipet.
3. Slowly pour the sample solution from step III.B.2, all at once, onto the resin bed.
4. Without touching the inside walls of the beaker, use the marked pipet to rinse down the walls of the sample beaker with 2 mL of 8 *M* HCl. Pour this rinse solution into the column. Repeat four more times for a total of five rinses.
5. Rinse the outside of the beaker spout into the column with a few drops of 8 *M* HCl.
6. Use the column effluent for appropriate analysis or discard if not needed.

D. Washing

1. Using the marked pipet, wash down the sides of the column reservoir with 5 mL of 8 *M* HCl. Allow the column to drain to the level of the resin bed.
2. Continue adding 8 *M* HNO₃ in 5 mL increments, allowing column to drain between additions, until 20 mL of wash have been added.

NOTE: If Cu and Co are present as impurities, Cu may not be separated completely with 8 *M* HNO₃ washes and Co will be separated only with 5 *M* HNO₃. If Cu and Co

are to be removed, wash the column with 15 mL of 8 *M* HNO₃, added in 5-mL increments, followed by 15 mL of 5 *M* HNO₃.

NOTE: Excess washing can result in loss of uranium.

3. Use the column effluent for appropriate analysis or discard if not needed.

E. Eluting

1. Place a clean, labeled 100-mL beaker under the column. The tip of the column should be approximately centered with respect to the beaker. If necessary, adjust the column holder to minimize distance between column tip and beaker rim.
2. Using the marked pipet, add 1 *M* HCl eluting solution to the column dropwise, washing the sides of the reservoir. Add no more than 5 mL at a time. Allow the liquid level to drain to the top of the resin.

NOTE: Addition of less than 5-mL volumes of the eluting solutions at a time to the reservoir can result in better exchange kinetics and therefore better recoveries.

3. Continue adding 1 *M* HCl eluting solution in increments of 5 mL, allowing the column to drain between increments, until 40 mL of the eluant have been added. Allow the column to drain until the solution stops dripping.
4. Rinse the tip of the column into the beaker with a few drops of eluting solution.
5. Discard the resin and column (see section VII.B).
6. Carefully place the beakers into an enamel-coated tray, move them and place on a hot plate.

F. Fuming

NOTE: The hot plate settings given below are usual values, but heating efficiency may vary. It is the analyst's responsibility to observe the heating process and adjust the settings as appropriate.

1. Set the hot plate to the "2" setting.
2. Place the sample beaker on the hot plate and allow the sample to dry down. Do not boil or bake the sample. This initial evaporation step should take about four hours. Remove the sample and allow it to cool.
3. Using a new disposable pipet, rinse down the walls of the beaker with 2 mL of concentrated HNO_3 .

NOTE: If the sample does not contain aluminum, also add 2 mL of 3 M H_2SO_4 .

4. Place sample beaker back on the hot plate and allow the sample to fume dry. This fuming step should take less than an hour if only HNO_3 is present. The fuming time will be longer if H_2SO_4 is used.

NOTE: Turning the heat up too soon or too high may cause the sample to "pop" with a resultant loss of sample of up to several tenths of a percent. Be aware that based on their location on the hot plate, individual samples may fume at different rates.

5. Cool the sample until it can be comfortably handled (approximately 5 minutes) by reducing the setting to "LOW" or "OFF". Remove samples from the hot plate.
6. Reset the hot plate to the "2" setting.
7. Using the same disposable pipet as in step III.F.3, again rinse down the walls of the beaker with 2 mL of concentrated HNO_3 .
8. Repeat steps III.F.4 and III.F.5.
9. After the sample cools to room temperature, cover the sample with Parafilm and store in a safe place.

NOTE: Do not stretch the Parafilm over the beaker since it may fail now or at a later time. Instead, fold the Parafilm over the rim and gently press it in place.

G. Documentation

Enter the date and the LIMS numbers of all the samples ion exchanged in the ion-exchange log notebook kept in the laboratory. Initial and date the entry. Record in the analyst's notebook that the ion exchange was performed; include procedure number and issue date, the resin lot used, and sample LIMS numbers. Note unusual conditions observed during ion exchange.

IV. CALCULATIONS

Not applicable

V. QUALITY ASSURANCE

To ensure quantitative recovery of uranium (and/or impurities) prior to analysis, known or blind standards, as appropriate, are ion exchanged concurrently with unknown samples. At least two knowns or blinds must be included with every set of unknown samples purified by this procedure. Analytical results may not be used if knowns or blinds are out of control. See applicable analysis procedure for details. The analyst should determine if the problem lies with the analysis or with the ion exchange; if the problem lies with the ion exchange, the set of samples must be discarded and the purification repeated with a new set of knowns/blinds and samples.

VI. TRAINING

A. Training Plan

The trainee observes a qualified analyst perform the ion exchange procedure. The trainee should practice the procedure until comfortable with the operations. The trainee must then demonstrate the procedure to a qualified analyst with sufficient technique to obtain quantitative recovery.

B. Qualification Plan

The trainee must perform the ion exchange procedure on sets of 6 known uranium samples on each of 2 days (total 12 samples). The ion exchange samples are analyzed by a qualified analyst using the appropriate procedure. Results (%RD) must be acceptable based on the capabilities of the analytical procedure.

C. Regualification Requirements

An analyst must purify at least one set (6 samples) of known uranium samples each year to maintain qualification. Results of analysis must be acceptable, as described in section VI.B. An analyst whose qualification has lapsed may requalify by meeting the same requirements.

VII. SAMPLE STORAGE AND DISPOSAL

A. Sample Storage

Always store the samples in a radioactive materials area.

B. Disposal

The effluent collected before loading the uranium fraction can be discarded into the sink with water. The effluent collected after sample loading and before eluting the uranium fraction, if not collected for impurity analysis, may contain toxic or radioactive elements. This effluent should be collected in 2-L (or smaller) polyethylene bottles and marked "LL-LRW" (low level liquid radioactive waste), listing all possible radioactive elements that may be present. These bottles are transferred into 55-gallon waste drums which will be disposed of as per NBL Standard Operating Procedures.

Collect the resin in a 2 L polyethylene container and dispose of it as Low-Level Radioactive Waste. All solid objects, including plastic pipets and ion exchange columns, should be disposed of in the suspect Solid Radioactive Waste (SRW, or "red-top") container placed inside the laboratory. Glass beakers should be cleaned and re-used. Broken glass should be disposed of into appropriately labeled containers placed in the lab.

Waste and scrap treatment is more fully described in the NBL ESH Manual and in NBL Standard Operating Procedures.³

VIII. REFERENCES

A. General References

1. U. I. Narayanan, P. M. Mason, J. P. Zebrowski, M. Rocca, I. W. Frank, M. M. Smith, K. D. Johnson, G. J. Orlowicz, and E. Dallmann, "The Quantitative Ion Exchange Separation of Uranium from Impurities," NBL-331 (March 1995).
2. ANL Parts List CMT-D1198 and Engineering Drawings CMT-C1198-3, -6, and -16; CMT-D1198-2, -4, and -5; and CMT-E1198-1.
3. NBL-SOP-ESH-04, "Waste Disposal Requisition Process."

B. Suggested Reading

1. "Dowex: Ion-Exchange," The Dow Chemical Company, Midland, Michigan, 1959.
2. H.O. Haug, "Gmelin Handbook of Inorganic Chemistry, Uranium Supplement, D3, Anion Exchange," Springer-Verlag, Berlin, 1982.
3. "Guide to Ion-Exchange," Bio-Rad, Catalog Number 140-9997.
4. "Sorbent Extraction Technology Handbook," N. Simpson, K.C. Van Horne, eds., Varian Sample Preparation Products, Harbor City, CA, 1993.
5. J. Korkisch, "Handbook of Ion Exchange Resins: Their Application to Inorganic Analytical Chemistry," Volume I and II, CRC Press, Boca Raton, FL, 1989.
6. I.M. Kolthoff and P.J. Elving, "Treatise on Analytical Chemistry," Part II, Vol. 9, .

APPENDIX B

STATISTICAL PLAN OF ANALYSIS

STATISTICAL PLAN OF ANALYSIS

Introduction

This project will test two ion exchange resins (Bio-Rad AG 1-X8 and MP-1) to determine which one will most closely achieve quantitative recovery of uranium. Two analysts will be involved in the project, each performing identical experiments, to test for analyst to analyst variation. Samples will consist of approximately 12 mg uranium. Glass and Teflon sample beakers will be used to determine the effect of sample beaker type on recovery. Half of the uranium samples will be contained in Teflon beakers and half will be in glass. The ion exchange method to be used will include an 8 M HCl load, an 8 M HCl wash and a 1 M HCl eluant. Eluted samples will be titrated using the low-level modified Davies and Gray procedure.

Uranium Samples

Uranium samples will be prepared from a single uranium solution made from New Brunswick Laboratory CRM 112-A. The solution will be aliquanted into both glass and Teflon 50 mL beakers. Note: the uranium solution aliquanted into the glass beakers (94NU0094-001-xxx) is from a different bottle (of the same solution) than that aliquanted into the Teflon beakers (94NU0094-008-xxx): the assumption is made that these solutions are identical. One analyst will aliquant all samples. Half of these samples will be spiked with a previously developed combination of impurities. All impurity spiking will be completed by one analyst. Samples will be identified as an numbered aliquant from either bottle 1 or bottle 8. For example, aliquant number 53 from bottle 1 will be assigned the identification number 001-053 (or just 1-53).

A total of seventy-four uranium samples will be aliquanted into glass beakers and seventy-two will be aliquanted into Teflon beakers. Some of these aliquants will be used for additional resin testing and direct titration without ion exchange. Thirty-six of the remaining aliquants will be randomly selected from each beaker type for ion exchange in the statistical plan. Half of the thirty-six in each type of beaker will be randomly selected to be spiked with the mixture of impurities. At this point four classifications of samples will be present, with 18 aliquants of each: uranium samples in glass beakers, uranium samples in Teflon beakers, uranium samples with impurities in glass beakers, and uranium samples with impurities in Teflon beakers. The individual samples in each of these four classifications will be randomly assigned to the two resins (nine samples of each classification to each resin). One sample from each of the nine will be designated as a spare; the remaining eight were randomly assigned to one of the two ion exchange analysts. Thus, each analyst will perform ion exchange on four samples from each of the four classifications described above.

Analyses

The following plan of analysis is for a single ion exchange resin. The plan will be performed as described for both resins by both analysts. Each set of ion exchanges will include both pure samples and those containing impurities. Pure and impure samples will be paired for statistical analysis within beaker type. Following ion exchange, each sample will be assayed by the low-level modified Davies and Gray titration. The number of titrations will have to be limited due to the availability of low-level blinds.

Two analysts will perform the ion exchanges, and a different analyst will perform all of the titration work for this project.

Each analyst assigned to ion exchange work will complete ion exchange of 10 samples each day for two days (a total of 20 samples). Of the 20 samples, 16 are designated for subsequent titration, two are designated spares, one is a blank and one a blank containing impurities. Thus, for a single resin, each analyst will complete eight sets (16 samples) of pure and impurity-added uranium samples; one sample designated as a spare on each day, and one blank on each day (either reagent blank or impurity-added blank).

ION-EXCHANGE RACK SET-UP

Order for ion exchanging

Day 1

U		
U		
U		
U		
U		
U		
U		
B		
U		
U _i		

ION-EXCHANGE RACK SET-UP CONT.

Day 2

U
U
U
U
U
B
U
U
U
U_i

U = pure uranium sample
U_i = uranium sample with added impurities
BI = reagent blank (empty beaker)
BI_i = blank with added impurities (impurities only, no U)

Of the 40 ion exchanged samples and blanks, 32 samples will be titrated for each resin, four are designated as spares and the four blanks, with and without impurities, will be saved for subsequent impurities measurements at a later time. The 32 samples, the four uranium-containing spares, and the ten non-ion-exchanged samples will be transferred to appropriate laboratory personnel for fuming and titration.

The titrations will be completed as listed in the titration section of this plan of analysis. As stated earlier, there are four sample spares for each resin (two with uranium only, two with uranium and impurities). This will allow for replacement of a sample designated for titration when the designated sample is determined to be inadmissible due to overtitration or other questionable event. Note that spares are designated for each solution-analyst, but not for beaker type. Titrated spares replacing a different beaker type will not be included in the analysis to detect beaker type differences.

Titration

Thirty-six samples processed with ion exchange and ten uranium samples that have not been ion exchanged will be fumed and titrated. Thirty-two ion exchanged samples and eight uranium samples are in the titration plan; the remaining samples are designated spares. Three days of titration will be designated for each resin.

Samples for titration will be blinded to the analyst (uranium only or uranium with Impurities, Ion exchanged or non-Ion exchanged). These will be titrated randomly but with pairs titrated on the same day (non-Ion exchanged samples are not paired). A listing of designated spares will be prepared for the analyst assigned to titrate these samples in order that a sample with like treatment replaces the lost sample.

Titration order:

Day 1:

K B SI S SI S SU SU SI B S SU SI SI SU SU B

Day 2:

K B SU SI S SU SI SU B S SI SI SI SU SU B

Day 3:

K B SU S SI SI SU SI SU B SI SU SU S SI SI B

K = Known*

B = Blind*

S = Sample with only uranium; non-Ion exchanged

SU = Sample with only uranium; Ion exchanged

SI = Sample with uranium and impurities; Ion exchanged

*Knowns and Blinds are quality control checks to ensure the Davies and Gray titration is performing adequately. These samples are not Ion exchanged.

APPENDIX C

STATISTICAL ANALYSIS OF URANIUM PURIFICATION PROJECT

REPORT OF DATA ANALYSIS

URANIUM SAMPLE PURIFICATION PROJECT

This project tested two ion exchange resins to determine which one would most closely achieve quantitative recovery of uranium, with and without the presence of impurities in the uranium solutions. Two ion exchange resins (Bio-Rad AG 1-X8 and MP-1) were tested. Ion exchange was performed by two analysts in order to test for analyst-to-analyst variation in performing the procedure. The ion exchange method used included an 8 M HCl load, 8 M HCl wash and a 1 M HCl elution.

Samples were aliquanted such that each sample contained approximately 12 mg of uranium. Samples from bottle 1 were aliquanted into glass beakers and those from bottle 8 were aliquanted into Teflon beakers. The use of Teflon beakers was included to test for variation in sample container type. The samples in Teflon beakers were aliquanted at a later date.

Following ion exchange, all samples were titrated using the low-level modified Davies and Gray procedure by one analyst with many years of experience in uranium titration. Spare samples were identified by resin, beaker type, and ion exchange analyst.

Uranium Samples

Uranium samples were prepared from a single uranium solution made from CRM 112-A. After preparation, this solution was stored in a series of glass bottles. The solution was aliquanted into appropriately-sized beakers (both glass and Teflon). One bottle of uranium solution (samples 94NU0094-001-xxx) was aliquanted into glass beakers before the Teflon beakers arrived from the supplier; the aliquants in the Teflon beakers were from another bottle (samples 94NU0094-008-xxx) of the same solution. The assumption was made that these bottles are the same. One analyst aliquanted all samples into the beakers, and spiked half of them with a previously prepared combination of impurities. Aliquants to be spiked were selected using a random number generator.

Seventy-four uranium samples were aliquanted into glass beakers and seventy-two were aliquanted into Teflon beakers. Some of these aliquants were used for additional resin testing, and some were subjected to direct titration without ion exchange. Thirty-six of the remaining aliquants were randomly selected from each beaker-type for ion exchange. Because half of the aliquants in each type of beaker were randomly selected to be spiked with a mixture of impurities, 18 aliquants from each of four classifications of samples were prepared: uranium samples in glass beakers, uranium samples in Teflon beakers, uranium samples

spiked with impurities in glass beakers and uranium samples with impurities in Teflon beakers. The samples in each classification were again randomly assigned to the two resins (nine samples of each classification to each resin). One sample from each of the nine was designated as a spare; the remaining eight aliquants of each type were randomly divided between the two ion exchange analysts.

In order to look at the effect of the impurities spike, spiked samples were paired with pure samples. Both aliquants of each numbered pair were taken from the same bottle of uranium solution, ion exchanged in the same ion exchange batch and titrated on the same day.

Ion Exchange and Titration

The ion exchange procedure was first completed for all AG 1-X8 samples by the two ion exchange analysts, each performing this procedure for two days. Designated pairs of aliquants were ion exchanged in the same batch, on the same day by the same analyst. Upon completion, the uranium samples and spares were transferred for low-level Davies and Gray titration. Additional samples from the same bottles used to aliquant the original solutions were added to the ion exchanged samples. These non-ion exchanged samples were randomly distributed in the three day plan for titration.

During preparation for titration, the ion exchanged samples for AG 1-X8 became contaminated and were deemed unusable. Ion exchange continued with samples using AG MP-1. The remaining aliquants taken from bottles 1 and 8 were used for the second round of ion exchange using AG 1-X8. Again, samples were randomly assigned to be spiked with the prepared solution of impurities. Because of the limited number of aliquants available, some of the samples originally designated to be titrated without ion exchange were replaced with aliquants from bottle 5, a bottle of the same solution used to prepare knowns for the low-level titration. The six non-ion exchanged samples to be titrated in each resin group consisted of three samples from bottle 1, two samples from bottle 5, and one sample from bottle 8.

Final ion exchange and titration procedures were completed on the dates listed in Table 1.

Table 1 Dates for Ion exchange and Titration Procedures		
	Ion Exchange	Titration
AG 1-X8	October 12-15, 1994	October 19-21, 1994
AG MP-1	October 4-7, 1994	October 11-13, 1994

Titration of Blind Samples

Low-level blinds were run at the beginning, middle, and end of each titration day according to the prescribed titration plan. Table 2 is a summary of the blinds data. All blinds were well within the 0.30 limit set for low-level Davies and Gray titrations. The percent relative differences (%RD) were positive for all blinds on the second titration day for AG 1-X8 and on the third day of titration with AG MP-1.

Table 2 Summary of Titration Blinds Data %RD is calculated from known uranium content.						
	AG 1-X8			AG MP-1		
Titration Day	n	Mean %RD	Std. Dev.	n	Mean %RD	Std. Dev.
Day 1	3	0.01040	0.08816	3	-0.02203	0.05470
Day 2	3	0.04043	0.03596	3	-0.00093	0.02699
Day 3	3	0.04127	0.04442	3	0.01050	0.01119

Knowns run at the beginning of each day were also well within the limits for low-level titration. These data were added to the non-ion exchanged data as they were prepared from the same solutions and subjected to the same analyses as the other non-ion exchanged samples. Thus the non-ion exchanged samples for each resin included three samples from bottle 1, five samples from bottle 5 and one sample from bottle 8.

Outlier Detection and Outliers in Samples

All outliers were in the negative direction for both resins, indicating a possible loss of sample upon ion exchange. Two were very obvious outliers with %RDs in excess of -1.2%, one from each resin. The outlier program detected four additional outliers from AG MP-1 and one additional outlier from AG 1-X8. Two outliers could be considered borderline (one from each resin) at -0.32% and -0.35%. The other three outliers detected by the outlier program in AG MP-1 were in the -0.6% RD range, all from same bottle (bottle 1, aliquanted into glass beakers), and all ion exchanged by the same analyst (one on ion exchange-Day 1 and two on ion exchange-Day 2). Thus five outliers were removed from the AG MP-1 data and two were removed from the AG 1-X8 data. Four of the outliers had been aliquanted into glass beakers and three were in Teflon beakers. Also, four of the outliers in AG MP-1 deleted two complete pairs of data.

Titration Data

After the removal of outliers and data without suitable spare replacements, 39 samples were left in the AG 1-X8 group (including 3 knowns and 6 other non-ion exchanged samples) and 36 samples were in the AG MP-1 group (including the same total number of non-ion exchanged samples). The overall results of these two groups are similar, as illustrated in Table 3.

Table 3 AG 1-X8 and AG MP-1 Titration Results NBL Modified Low-Level Davies and Gray Titrations						
%RD is calculated from known uranium content						
	AG 1-X8			AG MP-1		
	n	Mean %RD	Std.Dev.	n	Mean %RD	Std.Dev.
Non-ion exchanged	9	0.01184	0.04312	9	-0.03398	0.05169
Ion exchanged, Pure samples	16	-0.09191	0.10072	13	-0.02843	0.06799
Ion exchanged, Impurities added	14	-0.05081	0.06021	14	-0.06029	0.08057
All ion exchanged data	30	-0.07273	0.08548	27	-0.04495	0.07512
TOTAL	39	-0.05321	0.08527	36	-0.04220	0.06948

The difference between the results of the non-ion exchanged samples titrated with samples of each resin is in the same direction as the difference in the QC blinds shown in Table 2. Quality control blinds run with the AG MP-1 samples were slightly lower than those run with AG 1-X8, but not to the extent of the difference in the non-ion exchanged samples. Correction for differences in the quality control titrations affected only the fifth decimal places of the means of the non-ion exchanged samples.

No titration day-to-day variation was observed; no bottle-to-bottle variation was observed; no significant difference in results was observed between the glass and Teflon beaker preparations. However, it is noted that the mean of the data from the AG 1-X8 Teflon beakers is larger than that of the AG 1-X8 glass beakers when all samples are considered. However, when only the ion exchanged data are compared the means of the %RDs from titrations in both the glass and Teflon beakers are quite similar.

Ion Exchange Data

For each ion exchange resin, samples were ion exchanged by two analysts on two ion exchange days. Pairs of data were in the same ion exchange batch. Pairs were randomly assigned to analysts and batches. A summary of the %RDs from titration data, arranged by ion exchange analyst and ion exchange day, are illustrated in Table 4.

Table 4 Summary of Titration Data by Ion Exchange Analyst and Ion Exchange Day %RD is calculated from known uranium content						
	AG 1-X8			AG MP-1		
	n	Mean %RD	Std. Dev.	n	Mean %RD	Std. Dev.
ION EXCHANGE ANALYST 1						
IE-Day 1	7	0.00397	0.03382	7	-0.00933	0.07528
IE-Day 2	8	-0.04514	0.06804	7	-0.05270	0.07337
TOTAL	15	-0.02222	0.05872	14	-0.03101	0.07488
ION EXCHANGE ANALYST 2						
IE-Day 1	8	-0.08904	0.04378	7	-0.02995	0.04928
IE-Day 2	7	-0.16233	0.09453	6	-0.09495	0.08940
TOTAL	15	-0.12324	0.07887	13	-0.05995	0.07538
All ion exchange data	30	-0.07273	0.08548	27	-0.04495	0.07512

The ion exchange plan included eight samples (four pairs) to be ion exchanged on each day. The missing data represent the removed outliers. The following observations were made:

- In AG 1-X8, IE-analyst 2's overall mean %RD and mean %RDs for each ion exchange day are noticeably lower than those of IE analyst 1. The lowest AG 1-X8 %RD from IE analyst 1 data was -0.158 while the lowest for IE-analyst 2 was -0.260. All AG 1-X8 data for IE-analyst 2 have negative %RDs.
- The mean %RD of IE-analyst 2 data in AG MP-1 is also lower than that of IE-analyst 1, but the difference between means is not as great as that in the AG 1-X8 data. The spread of the data in AG MP-1 is similar for both analysts.

- c). The mean %RD for each IE-analyst in AG MP-1 is not significantly different from that of the non-ion exchanged data for AG MP-1. The mean %RD of the non-ion exchanged data for AG 1-X8 is positive though not significantly different from the data of IE-analyst 1 for AG 1-X8. However, the non-ion exchanged data for AG 1-X8 is significantly different from the IE-analyst 2 data for AG 1-X8.
- d). In AG 1-X8 IE-analyst 1 data, IE-day 1 has very precise data, while the IE-day 2 data is lower with considerably more spread (note that the standard deviation has more than double for a similar number of samples).
- e). The difference in ion exchange day data is also evident in AG MP-1 data--the mean %RD for IE-day 2 is consistently lower than the mean %RD for IE-day 1 (for both IE-analysts). However, in AG MP-1, IE-analyst 1 data have very similar standard deviations for both IE-days, while the AG MP-1 IE-analyst 2 data show that the standard deviation of IE-day 2 is considerably larger than that of IE-day 1. (Also note that five outliers were removed from AG MP-1 data.)

Significant Analysis of Variance Results

Analyst-to-analyst variations are significantly different for AG 1-X8 ($p = .0004$) but not for AG MP-1. Table 4 illustrates agreement between the analysts for IE-day 1 but not for IE-day 2. IE-day-to-IE-day variations (with analysts combined) are marginally significant for both resins ($p = 0.08$ and 0.06 , in that order).

Analysis of variance procedures, performed including both ion exchange analyst, ion exchange day and the interaction between IE-analyst and IE-day indicate significant contributions to the amount of variation in the %RD data for AG 1-X8 by both the ion exchange analyst ($p = 0.0002$) and the ion exchange day ($p = 0.015$), but only the ion exchange day has some influence (and only marginally significant with $p = 0.067$) on the amount of variation in AG MP-1. Interactions between IE-analyst and IE-day were not significant in either resin. These designated contributions to the variation in AG 1-X8 indicate that samples were ion exchanged with this resin by only one analyst on only one day may have somewhat biased results. Thirty-six percent of the variation in AG 1-X8 data may be explained by IE-analyst while less than four percent of the variation in AG MP-1 is explained by that variable. The model including both IE-analyst and IE-day shows that 50% of the variation in AG 1-X8 and 18% of the variation in AG MP-1 are explained by these two variables.

No titration day-to-day variation was observed; no bottle-to-bottle variation was observed; no significant difference was observed between the samples in glass and those in Teflon beakers.

Significant differences were observed among the "impurity types" of the AG 1-X8 samples. The three impurity types compared were the non-ion exchanged samples, the pure solution ion exchanged samples and the impurity spiked ion exchanged samples. The effect of this "impurities variable" was significant at the $p = 0.02$ level in AG 1-X8, but not at all significant in AG MP-1. When only the ion exchanged samples are considered, there are no observable significant differences between the uranium only and the spiked (with added impurities) samples in either resin. Thus the significant difference observed in AG 1-X8 is due to the results from the non-ion exchanged samples. In AG 1-X8, the group with the lowest %RDs is the uranium only solution ion exchanged samples group (mean %RD = -0.092), and the non-ion exchanged group has the highest %RDs (mean %RD = +0.118).

All three "impurity" groups have similar means in resin AG MP-1 data. Additional analyses of variance by impurity group indicate that resin AG 1-X8 and AG MP-1 are marginally different ($p = 0.06$) in the uranium only ion exchanged samples; the %RDs from AG MP-1 samples are generally higher (less negative) than those from AG 1-X8 samples. Both Resins 1 and 2 have similar results on the impurity-spiked ion exchanged samples. The non-ion exchanged samples are significantly different ($p = 0.03$), by resin group, as discussed earlier.

Paired Differences between Impurity-Spiked and Pure Solutions

A variable $DIFF = \%RD_{no\ imp} - \%RD_{imp\ added}$ was created for analyzing the matched pairs of samples from the same bottle. It is noted that the mean DIFF is negative (-0.041) for AG 1-X8 and positive (0.036) for AG MP-1. These two values were determined to be significantly different from each other with the Student t-test. Individually, the mean value of this difference is not significantly different from zero in the AG MP-1 data, but is marginally different from zero for AG 1-X8 ($p = 0.069$).

Comments and Selection of Ion Exchange Resin

Neither AG 1-X8 nor AG MP-1 clearly indicates a better recovery of uranium after ion exchange than the other resin. AG MP-1 produced more consistent results than those of AG 1-X8. However, AG MP-1 data had five outliers removed and the results of the non-ion exchanged samples titrated with resin AG MP-1 are still somewhat difficult to explain. The locations of the outliers removed from AG MP-1 data also raises additional questions.

The change in sign of the difference between the paired aliquants also should be considered. The uranium recovery of the samples containing impurities was much better than that of Uranium only ion exchanged samples with AG 1-X8. It would seem logical that it should be easier to recover the uranium from uranium

only samples than those containing uranium with impurities. The uranium only ion exchanged samples showed a better level of recovery (though not significantly better) after ion exchange on AG MP-1 than did the impurities-spiked samples with that resin.

Both resins produce biased results after titration: for AG 1-X8, the 95% confidence interval of the mean is 0.03192; thus the AG 1-X8 mean with confidence interval expression -0.07273 ± 0.03192 for all ion exchanged samples, illustrates the bias. The AG MP-1 mean with 95% confidence interval expression, -0.04495 ± 0.02972 for ion exchanged samples, illustrates a similar bias.

The titration results of samples prepared by ion exchange analyst 1 are similar for both resins. However, the titration results for the two resins were not similar for samples prepared by ion exchange analyst 2: the bias for AG 1-X8 is approximately double that of AG MP-1 for similar standard deviations. Titration results for samples from both ion exchange days of AG 1-X8 by analyst 2 indicate a significant negative bias. In all cases (both analysts, both resins), day 2 results have lower levels of recovery than those of IE-day 1.

Titration precision is similar for all titration days of AG MP-1 but not for AG 1-X8. The precision of titrations of samples prepared in both glass and Teflon beakers containing AG 1-X8 are similar, although those titrated without being ion exchanged, contained in glass beakers, exhibit less bias. Conversely, the bias of AG MP-1 samples in both types of beakers have similar means (bias), both with and without ion exchange, but the precisions are somewhat different depending on container type.

**SUMMARY TABLE OF %RD FOR URANIUM RECOVERY AFTER ION EXCHANGE FROM TWO RESINS
AND NON-ION EXCHANGED SAMPLES TITRATED CONCURRENTLY**

	AG 1-X8 Bio-Rad 1-X8			AG MP-1 Rio-Rad MP-1		
	n	Mean %RD	Std.Dev.	n	Mean %RD	Std. Dev.
ION EXCHANGED SAMPLES	30	-0.07273	0.08548	27	-0.04495	0.07512
IE-Day 1 Analyst 1	7	0.00397	0.03382	7	-0.00933	0.07528
IE-Day 2 Analyst 1	8	-0.04514	0.06804	7	-0.05270	0.07337
TOTAL IE-Analyst 1	15	-0.02222	0.05872	14	-0.03101	0.07488
IE-Day 1 Analyst 2	8	-0.08904	0.04378	7	-0.02995	0.04928
IE-Day 2 Analyst 2	7	-0.16233	0.09453	6	-0.09495	0.08940
TOTAL IE-Analyst 2	15	-0.12324	0.07887	13	-0.05995	0.07538
Non-ion-exchanged Samples	9	0.01184	0.04312	9	-0.03398	0.05169
Ion exchanged, Pure Samples	16	-0.09191	0.10072	13	-0.02843	0.06799
Ion exchanged, Impurities-added Samples	14	-0.05081	0.06021	14	-0.06029	0.08057
Titration Day 1	8	-0.11769	0.10703	6	-0.08305	0.08250
	12	-0.06899	0.11542	10	-0.06822	0.07421
Titration Day 2	10	-0.08235	0.08700	10	-0.04227	0.07724
	13	-0.05566	0.08339	13	-0.03937	0.07014
Titration Day 3	12	-0.04304	0.05784	11	-0.02660	0.06820
	14	-0.03741	0.05530	13	-0.02503	0.06429
GLASS	15	-0.07171	0.08580	14	-0.04453	0.05926
	23	-0.04071	0.08454	23	-0.04004	0.05665
TEFLON	15	-0.07368	0.08814	13	-0.04539	0.09059
	16	-0.07118	0.08574	15	-0.04940	0.08698
Paired Difference, $RD_{no. imp} - RD_{imp}$	14	-0.04121	0.07772	13	0.03647	0.08955

NOTE: The blocks for Titration Days and Glass/Teflon data include both ion exchanged and non ion exchanged data. The statistics for the smaller n in each grouping contains the ion exchanged only data.

DISTRIBUTION LIST

No. of Copies

Ames Laboratory
Ames, IA

ATTN: M. C. Edelson
Library

2

Argonne National Laboratory
Argonne, IL

ATTN: D. W. Green
R. B. Perry

2

Argonne National Laboratory
Idaho Falls, ID

ATTN: J. R. Krusul

1

Association Francaise de Normalisation
Paris, France

ATTN: J. C. Courtier

1

Babcock & Wilcox Company
Lynchburg, VA

ATTN: Information Services

1

Battelle Memorial Institute
Pacific Northwest Laboratory
Richland, WA

ATTN: Technical Information Section

1

Bendix Kansas City Division
Kansas City, MO

ATTN: L. Stratton

1

British Nuclear Fuels Ltd.
Capenhurst Works, United Kingdom

ATTN: Library

1

British Nuclear Fuels Ltd.
Springfield Works, United Kingdom

ATTN: Library

1

British Nuclear Fuels Ltd.
Windscale & Calder Works, United Kingdom

ATTN: Library
I. McGowan

2

Brookhaven National Laboratory
Upton, NY

ATTN: Research Library
Reference Section
J. Indusi
A. Reisman (ISPO)

4

Bundesanstalt für Materialprüfung Berlin, Federal Republic of Germany	ATTN: T. Tamberg	1
CEC-JRC Geel Institute for Reference Materials and Measurements Geel, Belgium	ATTN: P. DeBievre K. Meyer	2
CEC-JRC Karlsruhe Institute für Transurane Karlsruhe, Federal Republic Germany	ATTN: W. Bartscher	1
CEC-JRC Ispra Ispra, Italy	ATTN: M. Cuypers	1
COGEMA Center de la Hague Cherbourg, France	ATTN: Mme. Baude	1
COGEMA Etablissement de Pierrelatte Pierrelatte, France	ATTN: R. Picard	1
Comitato Nazionale per l'Energia Nucleare Rome, Italy	ATTN: A. Moccia V. Pagliai	2
Commissariat à l'Energie Atomique Grenoble, France	ATTN: L. Chaussabel J. Dettwiller	2
General Atomic Company San Diego, CA	ATTN: Library	1
Gesellschaft zur Wiederaufarbeitung von Kernbrennstoffen mbH Eggenstein-Leopoldshafen Federal Republic of Germany	ATTN: R. Berg	1
International Atomic Energy Agency Safeguards Analytical Laboratory Seibersdorf, Austria	ATTN: S. Deron G. Bagliano	2

International Atomic Energy Agency Vienna, Austria	ATTN: E. Kuhn	1
Kernforschungszentrum Karlsruhe GmbH Karlsruhe, Federal Republic of Germany	ATTN: Library	1
Lawrence Livermore National Laboratory Livermore, CA	ATTN: W. Ruhter	1
Los Alamos National Laboratory Los Alamos, NM	ATTN: J. Phillips J. Tape D. J. Temer Report Librarian	4
Martin Marietta Energy Systems, Inc. Oak Ridge National Laboratory Oak Ridge, TN	ATTN: J. Wade M. Poutsma	2
Martin Marietta Energy Systems, Inc. Y-12 Plant Oak Ridge, TN	ATTN: J. Dorsey	1
Monsanto Research Corporation Miamisburg, OH	ATTN: R. Goss R. Sellar Library	3
Westinghouse Materials Co. of Ohio Cincinnati, OH	ATTN: Library	1
Nuclear Fuel Services Erwin, TN	ATTN: Library	1
Nuclear Material Control Center Safeguards Analytical Laboratory Japan	ATTN: K. Iwamoto	2
Reaktor Brennelement Union GmbH Hanau, Federal Republic of Germany	ATTN: R. Saiger	1

Westinghouse Hanford Operations Richland, WA	ATTN: B. Baird H. Meznarich A. E. Schilling	3
Rockwell International Rocketdyne Division Canoga Park, CA	ATTN: Library	1
EG&G Rocky Flats Plant Golden, CO	ATTN: J. Weiss Library	2
Sandia National Laboratories Albuquerque, NM	ATTN: D. Miyoshi	1
Studiecentrum voor Kernenergie Centre d'Etude de l'Energie Nucleaire Mol, Belgium	ATTN: P. DeRegge	1
United Kingdom Atomic Energy Authority Atomic Energy Research Establishment Harwell, United Kingdom	ATTN: D. Crossley	1
United Kingdom Atomic Energy Authority Dounreay Nuclear Power Development Establishment United Kingdom	ATTN: J. L. Drummond	1
STRUGHOLD (TSK) Aero Medical Library USAF School of Aerospace Medicine Brooks AFB, TX	ATTN: Documents Librarian	1
U.S. Department of Energy Washington, D.C.	ATTN: DOE Library	1
U.S. Department of Energy Office of Safeguards and Security Washington, D.C.	ATTN: E. J. McCallum O. B. Johnson B. Todd D. W. Crawford G. D. Smith	5

U.S. Department of Energy
Albuquerque Operations Office
Albuquerque, NM

ATTN: L. Muller

1

U.S. Department of Energy
Chicago Operations Office
Argonne, IL

ATTN: T. Gradle
C. Pietri
J. Huagen
R. Lang

4

U.S. Department of Energy
Idaho Operations Office
Idaho Falls, ID

ATTN: B. Green

1

U.S. Department of Energy
New Brunswick Laboratory
Argonne, IL

ATTN: H. R. Martin
K. Lewis
R. D. Oldham
W. G. Mitchell
M. M. Lachman

5

U.S. Department of Energy
Oak Ridge, TN

ATTN: Technical Information Center

4

U.S. Department of Energy
Oak Ridge Operations Office
Oak Ridge, TN

ATTN: R. Green
G. Proco

2

U.S. Department of Energy
Pittsburgh Energy Technology Center
Pittsburgh, PA

ATTN: Technical Library

1

U.S. Department of Energy
Richland Operations Office
Richland, WA

ATTN: J. Spracklen

1

U.S. Department of Energy
Rocky Flats
Golden, CO

ATTN: B. Reece

1

U.S. Department of Energy Oakland Operations Office Oakland, CA	ATTN: J. Turner D. Wilhem	2
U.S. Department of Energy Savannah River Operations Office Aiken, SC	ATTN: T. Williams	1
U.S. Department of Energy Schnectady Naval Reactors Office Schnectady, NY	ATTN: B. M. Erickson	1
U.S. Environmental Protection Agency Las Vegas, NV	ATTN: Director	1
U.S. Nuclear Regulatory Commission Region II Atlanta, GA	ATTN: E. McAlpine	1
U.S. Nuclear Regulatory Commission Region III Glen Ellyn, IL	ATTN: L. Creed	1
U.S. Nuclear Regulatory Commission Region V Walnut Creek, CA	ATTN: J. Montgomery	1
Westinghouse Electric Corporation West Mifflin, PA	ATTN: Librarian	1
Westinghouse Hanford Company Richland, WA	ATTN: Document Control	1
Westinghouse Idaho Nuclear Company Idaho Falls, ID	ATTN: Technical Library	1
Westinghouse Savannah River Co. Savannah River Site Aiken, SC	ATTN: P. T. Deason, Jr. A. Gibbs M. K. Holland W. A. Spencer Technical Library	5